



**3M**

**Tegaderm™**

Wound Care Solutions

3M™ Tegaderm™ Silicone Foam Dressings

**Where impressively  
strong meets  
amazingly gentle.**

# Wound care starts with skin care.


Skin is the body's largest organ, making it one of the most vulnerable to forces like pressure, friction and shear. As a wound care clinician, your mission is to protect and maintain the integrity of each patient's skin – but it's a complicated job with many factors to consider.

Take wound care dressings. Many clinicians choose to use silicone foam dressings because they're more gentle to skin than standard dressings, decreasing the risk of Medical Adhesive-Related Skin Injury (MARSI).


## Pioneering the science of strong and gentle.

3M invented gentle-to-skin medical adhesives more than 50 years ago, and we continue to innovate solutions that provide consistent adhesion with easy removal to help minimize the risk of MARSI. Learn how to protect your patients from MARSI at [3M.com/MARSI](https://www.3m.com/MARSI).


In a colorimetric protein study that measured skin cell proteins left on an adhesive after removing the adhesive from the skin, 3M's silicone adhesive removed **significantly fewer skin cells** than a standard acrylate adhesive – helping to reduce the risk of skin trauma.<sup>1</sup>

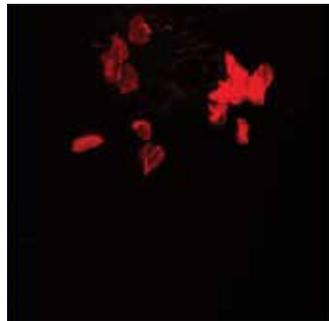

Skin cell proteins on 3M's gentle silicone adhesive.

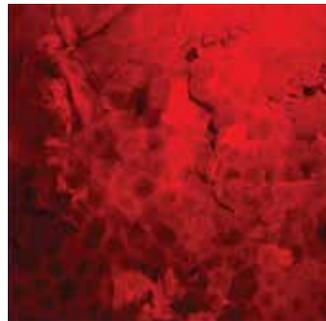

Skin cell proteins on a standard acrylate adhesive.

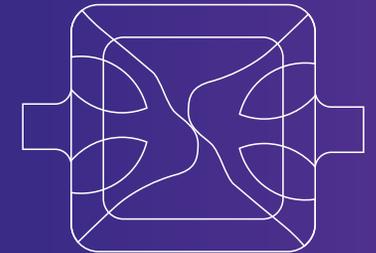


Not all silicone foam dressings are created equal. Unlike many competitors who license their adhesive technology, 3M scientists developed our unique silicone adhesive in our labs – applying their decades of adhesive expertise and innovation to your toughest dressing challenges.


3M™ Tegaderm™  
Silicone Foam  
Dressings feature  
layers of innovation.

**Patented spoke  
delivery system**  
enables easy application,  
even with gloves on, so  
your other hand is free  
to position the patient.

**Exclusive  
3M™ Tegaderm™  
film backing**  
combines 3M's  
adhesive innovation  
and film expertise  
in a breathable  
dressing cover.

**Moisture  
control layer**  
helps maintain  
an optimal  
moisture balance  
by facilitating  
evaporation through  
the film backing.

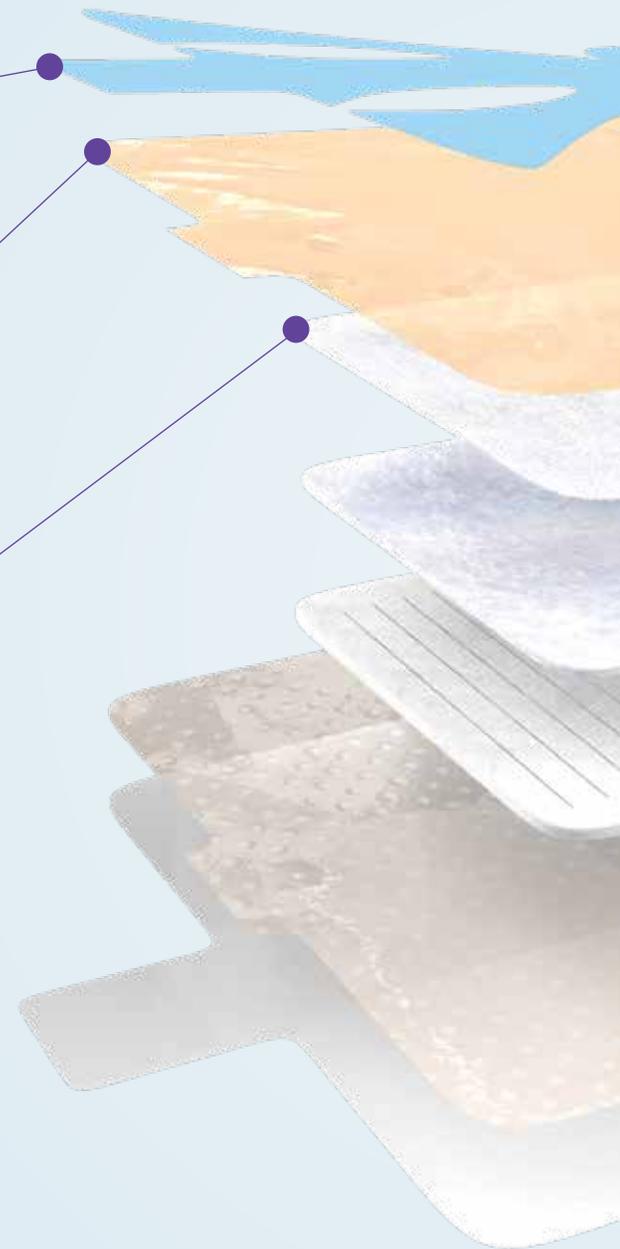

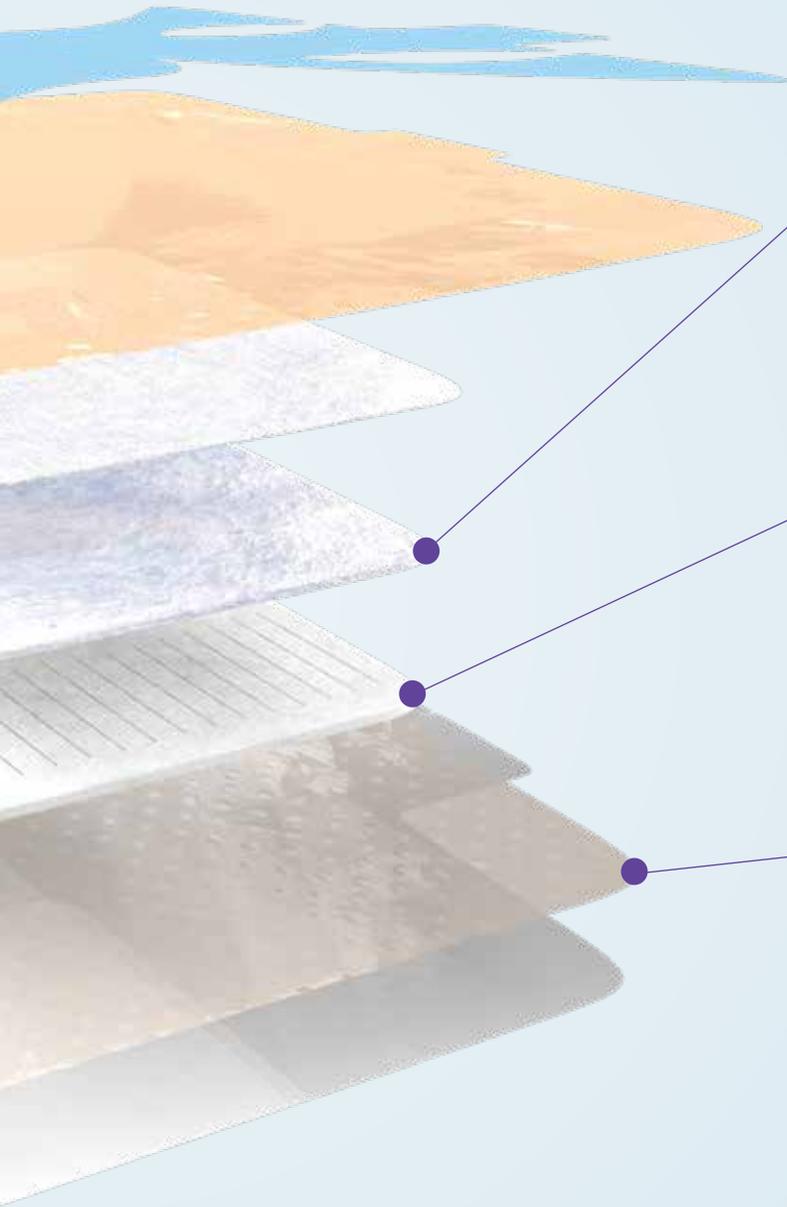

**Superabsorbent layer**  
further pulls moisture away  
and locks it in place to help  
minimize the backward  
migration that can cause  
periwound maceration.

**Fenestrated foam layer**  
improves flexibility and  
moisture absorption to help  
minimize the disruption of  
healing tissue.

**Proprietary  
silicone adhesive**  
keeps the dressing in place  
without damaging the  
wound area while allowing  
exudate to pass through.



## Your challenge:

Silicone foam dressings that don't wear as long as you need.



3M solution:

# Significantly longer wear time<sup>2</sup> plus gentle adhesion

3M™ Tegaderm™ Silicone Foam Dressings offer significantly longer wear time than the leading competitive silicone foam dressing<sup>2</sup> while being gentle to the skin. Which may help save your facility time and money on unscheduled dressing changes.

**2X<sup>3</sup>**  
longer  
wear time



**3M™ Tegaderm™  
Silicone Foam Dressing**



**Mepilex® Border Foam  
Dressing**

The 3M™ Tegaderm™ Silicone Foam Dressing wore 2X longer than the leading competitive silicone foam dressing when worn for seven days and lifted daily.<sup>3</sup>

2. 3M data on file. 4x4 and 6x6 dressings, based on *in vivo* studies EM-13977 and EM-13978.

3. 3M data on file. Median wear times. 6x6 dressings, based on *in vivo* study EM-13978.



## Your challenge:

Silicone foam dressings that lift and roll up in high-shear locations.



3M solution:


## Thin, low-profile edge

3M™ Tegaderm™ Silicone Foam Dressings are highly conformable and feature a thin, low-profile edge, helping to minimize the rolling and lifting that can impact adhesion and wear time.

## Dressing edge comparison



3M™ Tegaderm™ Silicone Foam Dressing

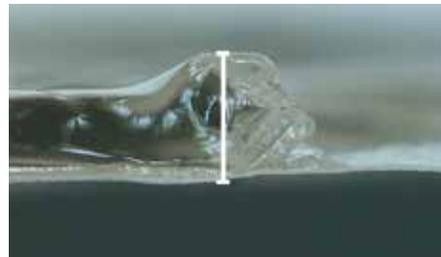

Mepilex® Border Foam Dressing

## Your challenge:

Wound drainage that pools under the dressing, causing periwound maceration.




3M solution:

## Unique multi-layer design

3M™ Tegaderm™ Silicone Foam Dressings have a unique multi-layer design that absorbs and evaporates moisture to help reduce the potential for skin maceration.

## Dressing saturation comparison


In a simulated invitro study with a highly exuding wound model under compression, the 3M™ Tegaderm™ Silicone Foam Non-Bordered Dressing had **significantly less fluid pooling** on the wound side of the dressing compared to the leading silicone foam dressing competitor.<sup>5</sup>

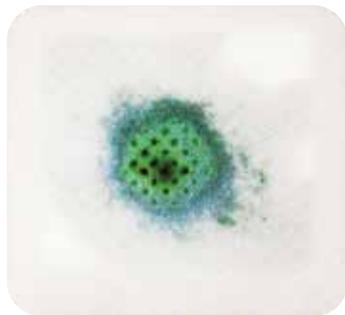

3M™ Tegaderm™ Silicone Foam Non-Bordered Dressing



Mepilex® Foam Dressing

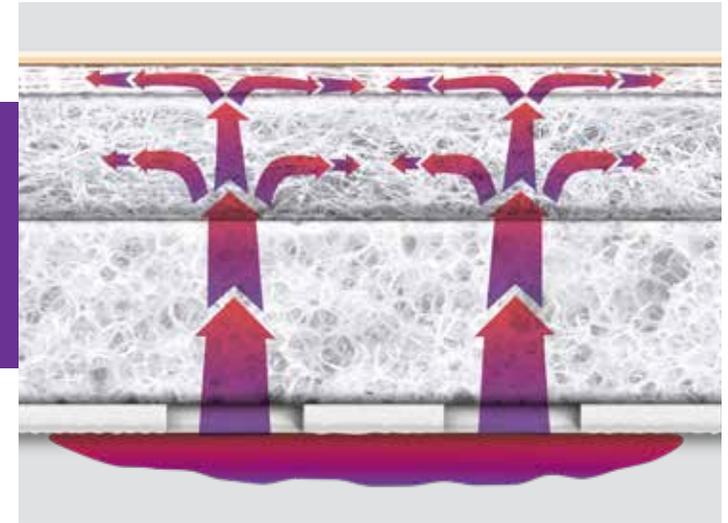

## Fluid management simulation

Our innovative layer technology is designed with a superabsorber that helps minimize the backward moisture migration that can cause maceration.

## Your challenge:

Anatomical locations that make dressing application difficult.



3M solution:

## Easy application

Our patented spoke delivery system enables easy application in challenging locations like the sacrum, for a more positive clinician experience.



Most evaluators (31 out of 33) were **“very satisfied”** or **“satisfied”** with ease of application for the 3M™ Tegaderm™ Silicone Foam Sacral Dressing.<sup>6</sup>

See how 3M™ Tegaderm™ Silicone Foam Dressings can make a difference in your practice. Visit [3M.com/SiliconeFoam](https://www.3m.com/SiliconeFoam) to request a free sample.



# Where smart protection meets pressure ulcer/injury prevention.

Facility-acquired pressure ulcers/injuries are a growing healthcare problem.<sup>7</sup> Not only can they lead to longer hospital stays and higher rates of readmission, but they can contribute to greater patient pain and suffering – and in some cases, premature mortality.<sup>8</sup>

As part of a comprehensive pressure ulcer/injury prevention plan, the use of polyurethane foam dressings to protect bony prominences from friction and shear should be considered to decrease the risk of pressure ulcer/injury development.<sup>9</sup>



7. Zaratkiewicz, S., Whitney, J. D., Lowe, J. R., Taylor, S., O'Donnell, F., & Minton-Foltz, P. (2010). Development and Implementation of a Hospital-Acquired Pressure Ulcer Incidence Tracking System and Algorithm. *Journal for Healthcare Quality*, 32(6), 44-51.

8. Health Research & Educational Trust (2016, January). Hospital Acquired Pressure Ulcers (HAPU) Change Package: 2016 Update. Chicago, IL: Health Research & Educational Trust. Accessed at [www.hret-hen.org](http://www.hret-hen.org).

9. European Pressure Ulcer Advisory Panel, National Pressure Injury Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers/Injuries: Clinical Practice Guideline. Emily Haesler (Ed.). EPUAP/NPIAP/PPPIA: 2019

## Ideal properties for dressings used in pressure ulcer/injury prevention:

When selecting a dressing for pressure ulcer/prevention, there are several ideal properties to consider – including the following from the National Pressure Ulcer Advisory Panel (NPUAP) *Prevention and Treatment of Pressure Ulcers: Clinical Practice Guidelines* document:

- ✓ **Ability to access and assess skin**  
Body areas at high risk for pressure ulcer/injury should be inspected often to detect early signs of pressure damage.

*Look for a dressing that can be lifted and re-adhered frequently for assessment without damaging the skin.*

- ✓ **Ease of application and removal**  
Applying and removing dressings from locations such as the heel and sacrum can be challenging, often requiring assistance to properly position the patient.


*Look for a dressing designed to make application easier, which in turn can help lead to fewer dressing failures and fewer unnecessary dressing changes.*

- ✓ **Ability to manage microclimate**  
Warm, moist skin is more vulnerable to the damaging effects of pressure and shear, which are recognized risk factors for pressure ulcer/injury formation.<sup>10</sup>

*Look for a dressing with properties that reduce the amount of moisture trapped at the skin's surface.*



- ✓ **Correct dressing size for high-risk locations**  
Anatomical sites that overlay a bony prominence, such as the heel and sacrum, account for more than half of all pressure ulcers/injuries<sup>9</sup> due to their vulnerability to pressure, friction and shear.

*Look for a dressing that is specifically designed for these high-risk locations and available in sizes to accommodate a wide range of body types.*

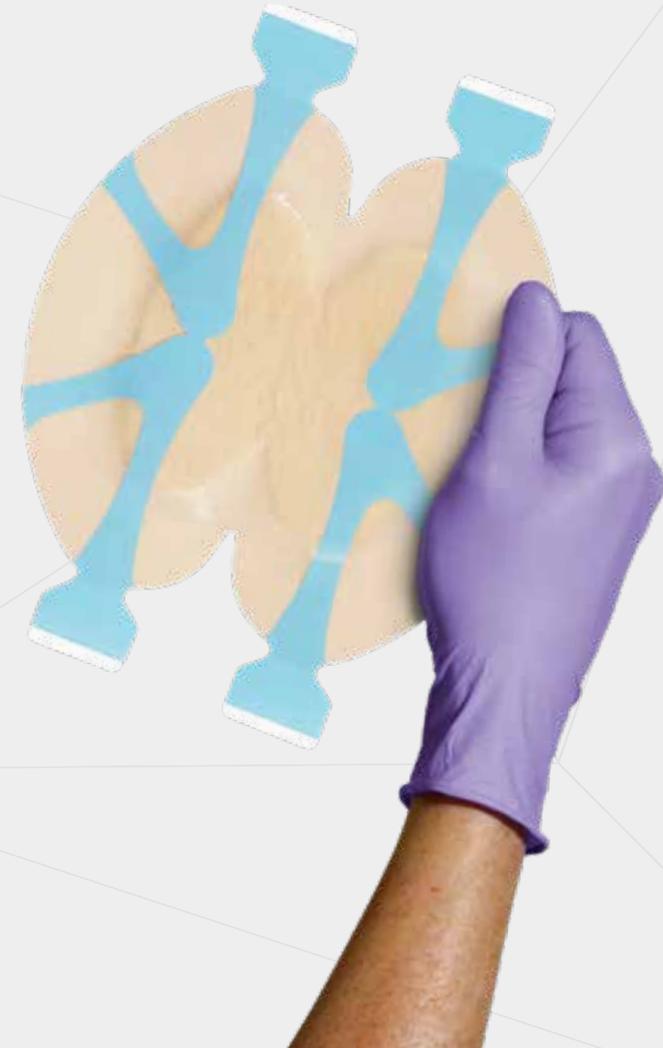


# Where variety meets value.

Choose 3M™ Tegaderm™ Silicone Foam Dressings as part of your wound management and pressure ulcer/injury prevention programs.

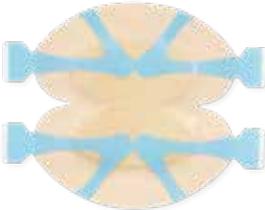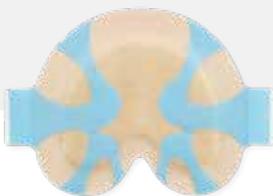	Product	Cat Number	HCPCS Code	Size	Dressings/Box	Boxes/Case
	Non-Bordered Dressing	90631	A6210	4 in. x 4.25 in. 10 cm x 10.8 cm	10	4
	Non-Bordered Dressing	90632	A6210	6 in. x 6 in. 15 cm x 15 cm	10	4
	Bordered Dressing	90643	A6212	2 in. x 2 in. 5 cm x 5 cm	10	6
	Bordered Dressing	90640	A6212	3 in. x 3 in. 7.5 cm x 7.5 cm	10	6
	Bordered Dressing	90641	A6212	4 in. x 4 in. 10 cm x 10 cm	10	6
	Bordered Dressing	90642	A6213	6 in. x 6 in. 15 cm x 15 cm	10	4
	Heel & Contour	90646	A6212	6.5 in. x 6.5 in. 16.5 cm x 16.5 cm	5	4
	Small Sacral	90647	A6213	6 in. x 6.75 in. 15 cm x 17 cm	10	4
	Large Sacral	90648	A6213	7.25 in. x 8.75 in. 18.5 cm x 22 cm	5	4

Visit [3M.com/SiliconeFoam](https://www.3M.com/SiliconeFoam) to learn more and request a free product sample.

HCPCS codes have been provided to assist you in the preparation of insurance claims. Please note, however, that the reimbursement information provided by 3M Health Care and its representatives is intended to provide general information relevant to coverage and coding for 3M products. Insurers reimbursement policies can vary and the use of the codes discussed here does not guarantee that an insurer will cover or pay at any particular level. Health care providers should exercise independent clinical judgment in choosing the codes which most accurately describe the products provided.

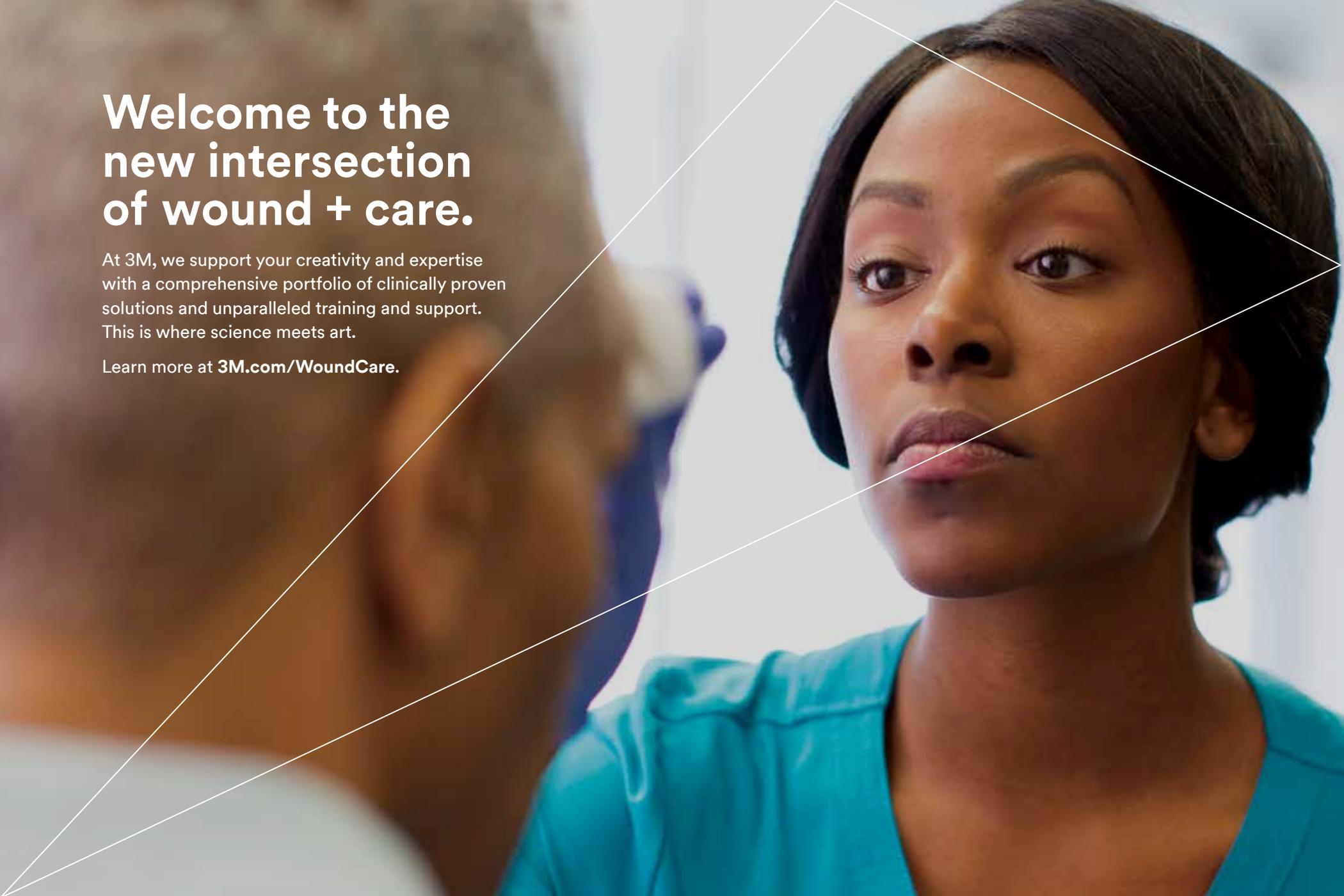
A woman with dark hair, wearing a teal scrub top, is looking intently at a patient whose hand is visible in the foreground. The background is a blurred hospital setting. A white geometric shape, resembling a triangle, is overlaid on the right side of the image.

# Welcome to the new intersection of wound + care.

At 3M, we support your creativity and expertise with a comprehensive portfolio of clinically proven solutions and unparalleled training and support. This is where science meets art.

Learn more at [3M.com/WoundCare](https://www.3M.com/WoundCare).

Discover why 3M™ Tegaderm™ Silicone Foam Dressings are an excellent choice for your wound management and pressure ulcer/injury prevention programs. Visit [3M.com/SiliconeFoam](https://www.3m.com/SiliconeFoam) to request a free sample.

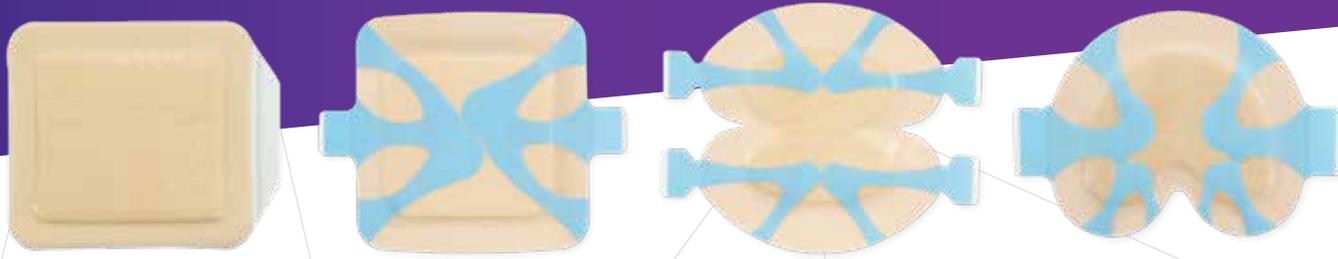

**3M**

**3M Company**  
2510 Conway Ave.  
St. Paul, MN 55144 USA

Phone 1-800-228-3957  
Web [3M.com/medical](https://www.3m.com/medical)

© 3M 2021. All Rights reserved. 3M and the other marks shown are marks and/or registered marks. Unauthorized use prohibited. 3M marks used under license in Canada. All other marks are property of their respective owners. 70-2011-6792-4