



3M Science.
Applied to Life.™

3M Biopharmaceutical Clarification and Purification for mAb/Recombinant Protein

Innovation for streamlined and accelerated workflows.

https://www.3mindia.in/3M/en_IN/bioprocessing-in/





One step can change everything.

A single step in the right direction can change more than just your perspective on bioprocessing — it can change your entire perception of what's possible. One step can transform your recombinant protein manufacturing process — from adopting 3M's next gen, single-stage products, to upgrading your process for greater scalability, efficiencies and faster speed to market.

Backed by decades of expertise in filtration and separation technologies, 3M's clarification and purification portfolio for mAb/recombinant protein solutions is designed to condense steps in the process while helping to improve efficiencies.

Reimagine bioprocessing

Standard bioprocessing lines can be fraught with challenges and outdated technology that can't keep up with new developments in higher-titer cell culture. From greater impurities (contaminants) and scale-up issues to prolonged downtime for cleaning and downstream bottleneck, we have a portfolio of solutions that help your line run smoother and faster with excellent purity and yields.

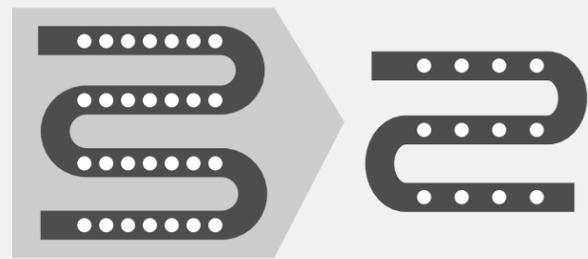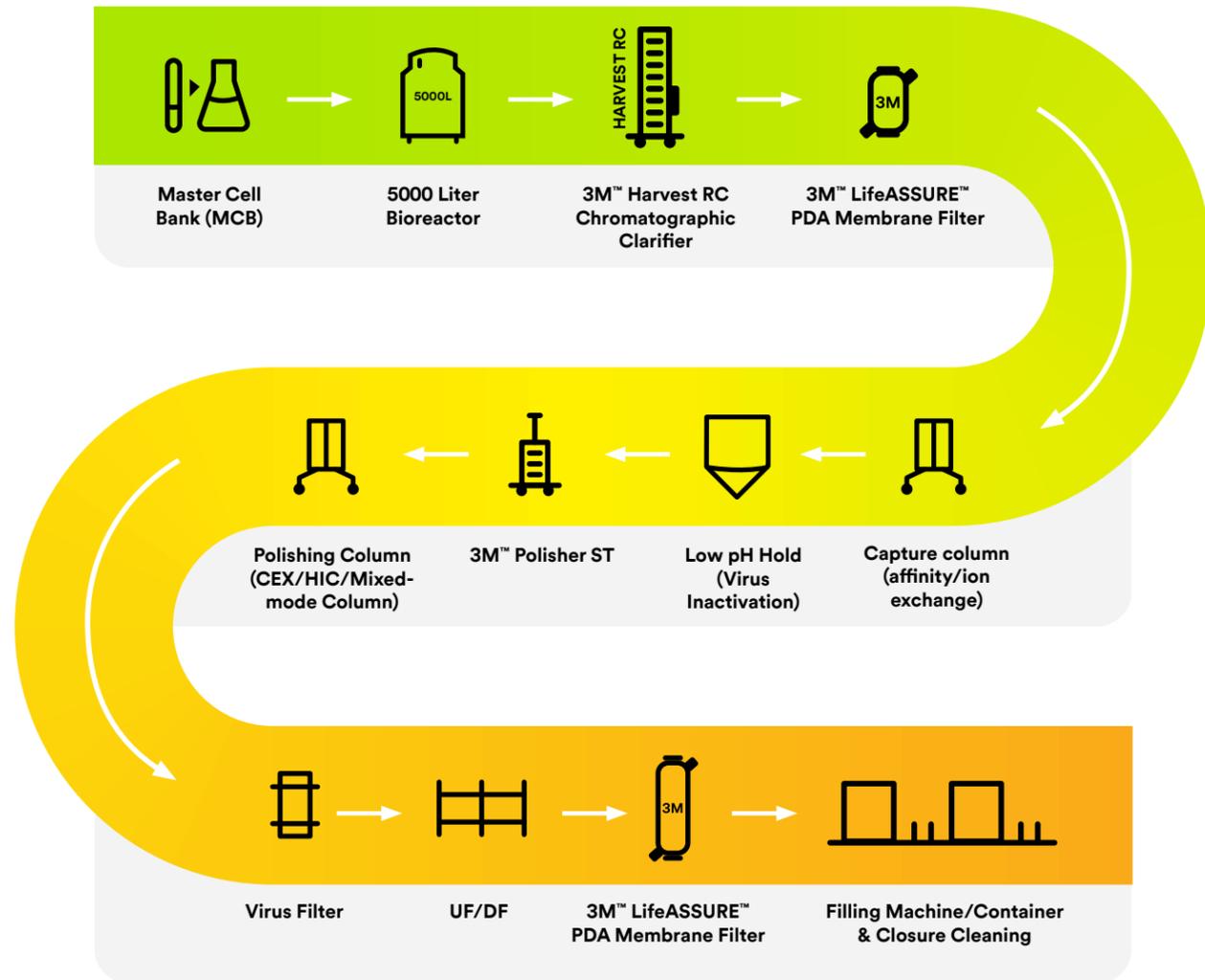

3M's portfolio for mAb/recombinant protein clarification and purification can take a standard bioprocessing line and condense several steps into one, saving you valuable time without sacrificing quality.

Advanced harvest, clarification and purification technologies

From lab scale-up to commercial manufacturing, we work with you across the entire process train for harvesting and purification to meet the ever-increasing demand to deliver high-quality biopharmaceutical products at an accelerated speed to market.

Harvest and Clarification Solutions



We offer a comprehensive range of harvest and clarification technologies — from standard filtration to next-gen chromatographic clarification — to help you achieve high yield and purity.

Products:

- 3M[™] Zeta Plus[™] Depth Filter
- 3M[™] Harvest RC Chromatographic Clarifier
- 3M[™] Emphaze[™] AEX Hybrid Purifier

Purification and Contaminant Removal Solutions



Our high-performance polishing and sterile filtration products are ideal for scale-up — delivering advanced purification and contaminant removal. We offer a single-use chromatographic polisher with functionalized membrane technology for efficiency and speed.


Products:

- 3M[™] Polisher ST
- 3M[™] LifeASSURE[™] PDA Membrane Filter



Process compression for intensified mAb manufacturing

Increasingly, you're expected to do more with your manufacturing process – including delivering with greater speed, within budget and under ever-growing challenges. 3M's single-use technology helps you maximize your process train's potential with our high performance single-stage clarification and purification solutions that include intensification options.

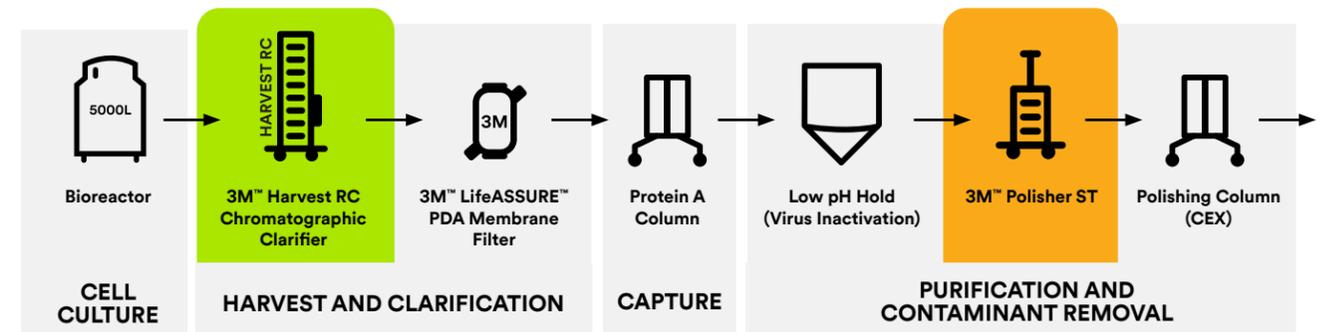


APPLICATION	FEATURED PRODUCTS
<ul style="list-style-type: none"> Harvest and Clarification Purification and Contaminant Removal 	<ul style="list-style-type: none"> 3M™ Harvest RC Chromatographic Clarifier Download product brochure 3M™ Polisher ST Download product brochure
TECHNOLOGY	KEY BENEFITS
Single-use fiber chromatography	Process compression, high yield, high purity

ADDITIONAL BENEFITS	
<ul style="list-style-type: none"> Streamlined handling of higher-concentration cultures (between 5-8% packed cell volume) Efficient separation of cell, cell debris and DNA from harvest fluid Consistent high yield Efficient removal of HCP 	<ul style="list-style-type: none"> Process compression – reduced unit operations Improved total cost of manufacturing compared to traditional setups 3M™ Polisher ST replaces the traditional AEX column 3M™ Polisher ST can replace multiple process steps

Optimizing opportunities

Whether you're building a new bioprocessing line or streamlining an existing one, we apply our 30 years of purification expertise to help optimize your operation with system-compatible products. Our innovative technologies are designed to improve upon standard bioprocessing by accelerating and condensing steps, leading to increased purity and yield.

See the 3M difference:

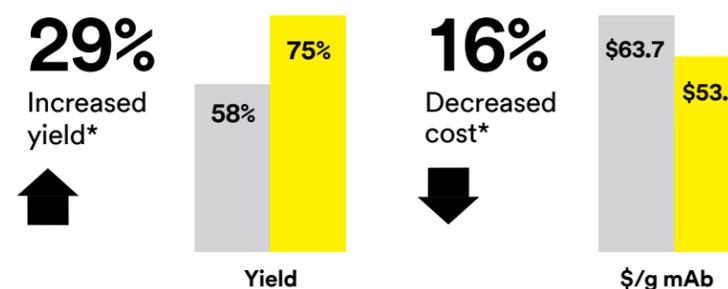


■ Legacy process with 2 stages of depth filters in harvesting & traditional AEX & CEX columns in polishing
■ 3M™ Harvest RC Chromatographic Clarifier and 3M™ Polisher ST

*Commercial Biosolve Process™ software package from Biopharm Services limited was used for modeling; The results are generated for the hypothetical scenarios. Depending on product and process conditions, the effects shown here may not be realized or may not be realized to the degree shown in the scenarios. 3M data on file.

Enhanced clarification for lines with centrifuge


3M's robust, single-use clarification solutions can help you improve purity while increasing efficiency on your traditional processing line. Our reliable 3M™ LifeASSURE™ PDA Membrane Filters help your control bioburden with fast flow rates. And our depth filters provide you with a full range of pore sizes and media types while eliminating the need for cleaning and validation steps — shortening setup in encapsulated systems.

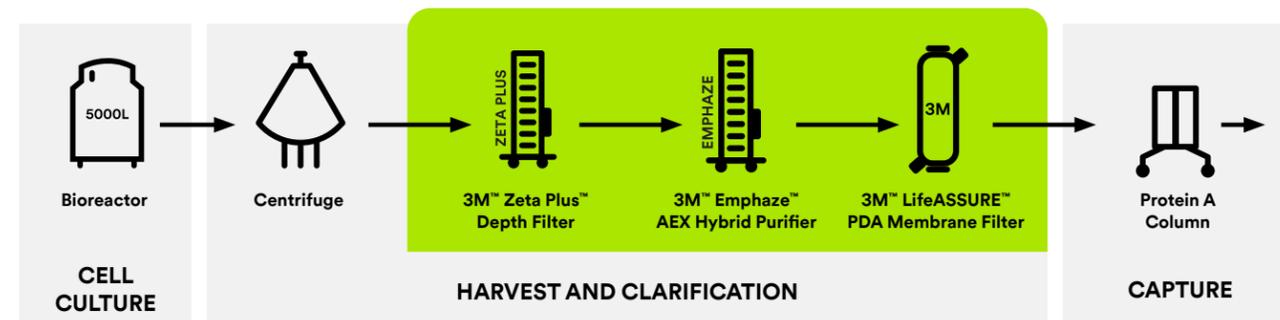

APPLICATION	FEATURED PRODUCTS
Harvest and Clarification	<ul style="list-style-type: none"> 3M™ Zeta Plus™ Depth Filter Download product brochure 3M™ LifeASSURE™ PDA Membrane Filter
TECHNOLOGY	KEY BENEFITS
Absorptive depth filter, micro porous membrane technology	HCP and turbidity reduction pre and post Protein A column
BENEFITS OF 3M™ ZETA PLUS™ DEPTH FILTERS	BENEFITS OF 3M™ LIFEASSURE™ PDA MEMBRANE FILTERS
<ul style="list-style-type: none"> 3 types of media: SP, LA, ZB Dual layer media enhances contaminant holding capacity ZB range has robust feed quality with higher charge level than SP and LA media 	<ul style="list-style-type: none"> High contaminant capacity Gamma compatible and sterile grades

Increased efficiency for lines with no centrifuge


Help improve your DNA and HCP reduction and increase product purity while reducing cell shear when you add one or two of our depth filters (depending on bioreactor size) and by adding our single-use chromatographic products. You may even be able to skip the second depth filter depending on your concentration of particles — saving steps and increasing process performance.

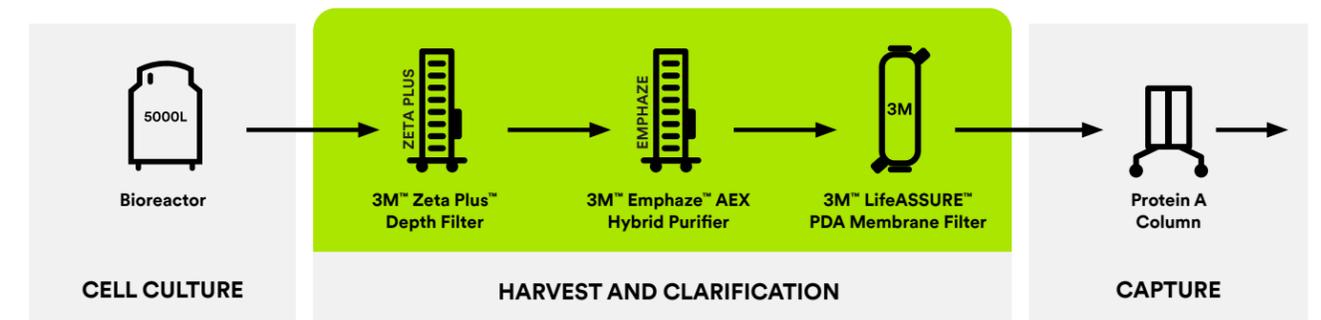

APPLICATION	FEATURED PRODUCT
Harvest and Clarification	3M™ Emphaze™ AEX Hybrid Purifier Download product brochure
TECHNOLOGY	KEY BENEFITS
Single-use fiber chromatography	Reduced cell shear and process variability
BENEFITS OF 3M™ EMPHAZE™ AEX HYBRID PURIFIER	
<ul style="list-style-type: none"> >10x less HCP post-Protein A >1,000x DNA reduction post-Protein A >10x less residual contaminants on Protein A column Improve clarified cell culture stability 	<ul style="list-style-type: none"> Reduce challenges for polishing train Improved process variability by eliminating centrifuge Reduce cell shear No need to clean and validate centrifuge

Process compression with single-stage harvest and clarification


You can eliminate the centrifuge in your line and condense steps when you utilize it with 3M's single-stage chromatographic products for harvest and clarification. It simultaneously provides whole cell removal and DNA reduction for high purity.

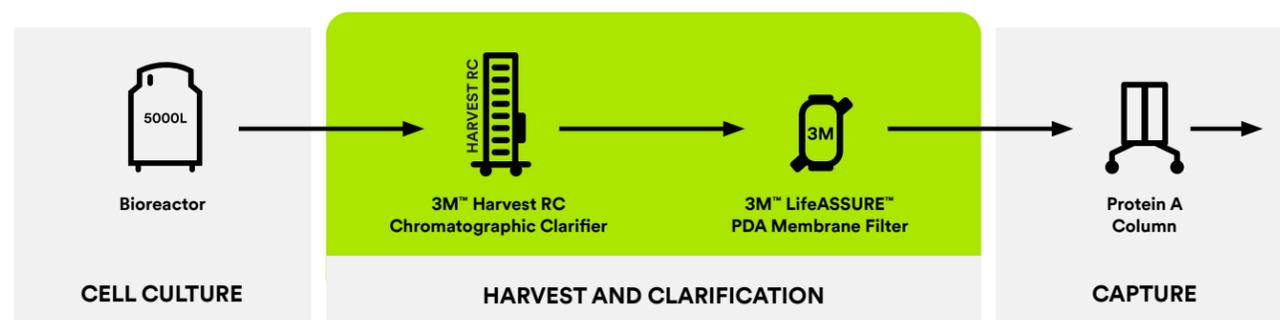

APPLICATION	FEATURED PRODUCT
Harvest and Clarification	3M™ Harvest RC Chromatographic Clarifier Download product brochure
TECHNOLOGY	KEY BENEFITS
Single-use fiber chromatography	Process compression, turbidity reduction, DNA removal
BENEFITS OF 3M™ HARVEST RC CHROMATOGRAPHIC CLARIFIER	
<ul style="list-style-type: none"> >10x less HCP post-Protein A High mAb product recovery (>95%) Turbidity reduction (<15NTU) DNA reduction (<500ppb) 	<ul style="list-style-type: none"> Process intensification Consistent high yield Predictable scalability

Process compression for polishing stage

Under certain operating conditions, you can compress multiple polishing steps into one with our chromatographic polisher — without sacrificing fluid purity. Our high capacity, yet cost-effective polishing solution replaces the traditional AEX column — providing robust viral clearance and DNA and HCP removal to help you overcome scale limitations.

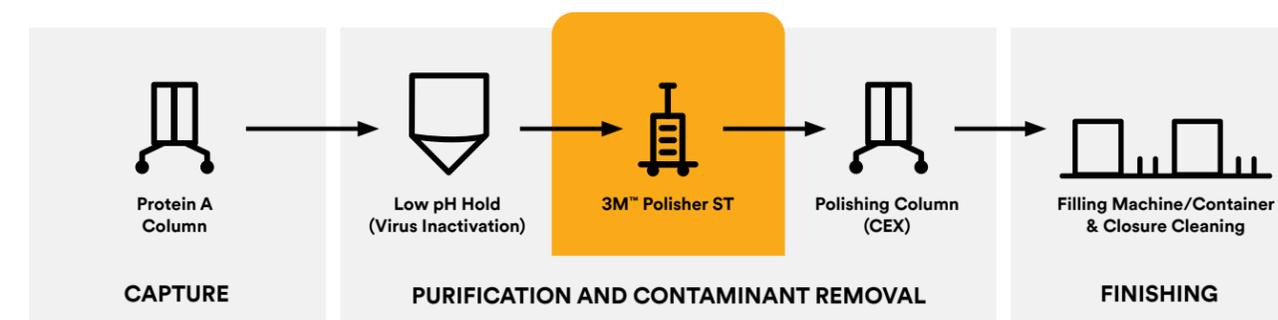

APPLICATION	FEATURED PRODUCT
Purification and Contaminant Removal	3M™ Polisher ST Download product brochure
TECHNOLOGY	KEY BENEFITS
Single-use fiber chromatography with functionalized membrane technology	Versatility, viral clearance, problematic HCP reduction, process compression
BENEFITS OF 3M™ POLISHER ST OVER TRADITIONAL AEX COLUMNS	
<ul style="list-style-type: none"> HCP reduction >50% above pH 5 Virus reduction >4 logs DNA clearance down to less than 10 ppb, or reduction >4 logs or below detection Platformability — works in a wide range of operating conditions 	<ul style="list-style-type: none"> High salt and low pH tolerance Lower buffer consumption Enables design of a more condensed polishing train Decreased total cost of ownership



Step into the future with our advanced fiber chromatographic technology

From lab to production scale, the innovative 3M™ Emphaze™ AEX Hybrid Purifier, 3M™ Harvest RC Chromatographic Clarifier and 3M™ Polisher ST help you streamline your recombinant protein production process, so you can deliver increased batch yields with predictable quality — saving time and money.

PRODUCT	3M™ Harvest RC Chromatographic Clarifier	3M™ Emphaze™ AEX Hybrid Purifier	3M™ Polisher ST
TECHNOLOGY PLATFORMS	<ul style="list-style-type: none"> Quaternary amine (Q) functionalized polypropylene nonwoven 0.2 μm polyether sulfone (PES) membrane 	<ul style="list-style-type: none"> Quaternary amine (Q) functionalized polypropylene nonwoven Highly asymmetric 0.2 μm polyamide membrane 	<ul style="list-style-type: none"> Quaternary amine (Q) functionalized polypropylene nonwoven Guanidinium (Gu) functionalized polyamide membrane
MECHANISM	AEX		
MATERIALS	Synthetic		
FEED STREAM	Direct harvest	Post depth filter, post centrifuge	Post Protein A, post-viral inactivation
SCALABILITY	Linear scalability from lab to multi-capsule commercial production		
APPLICATIONS	<ul style="list-style-type: none"> Direct harvest in a single stage Replacement of primary and 2nd depth filters and guard filters Replacement of centrifuge 	<ul style="list-style-type: none"> Replacement of 2nd depth filters in harvest and clarification Guard filter to Protein A Bioburden control 	<ul style="list-style-type: none"> Replacement of AEX column Replacement of depth filter and membrane filter prior to AEX column
FEATURES	<ul style="list-style-type: none"> High mAb product recovery (Capsules: >95%; Conical Tube and Well Plate >90%) Consistent cell loading capacity Turbidity reduction (<15NTU) DNA reduction (<500 ppb) Cell shear minimization 0.1 & 0.2 μm sterile filter protection 	<ul style="list-style-type: none"> Nominal 20-40% HCP and greater than 4 logs DNA reduction Consistent output turbidity (<5NTU) Typically, ≥ 4 LRV endotoxin reduction at low conductivity and neutral pH depending on buffer 	<ul style="list-style-type: none"> Operates in pH 5.0 to 9.0 and conductivity 3 to 20 mS/cm process conditions HCP reduction >50% above pH 5 Virus reduction >4 logs DNA clearance down to less than 10 ppb or greater than 4 logs or below detection Post-use integrity testable

Portfolio Overview

One step can improve your bioprocessing line.



3M™ Zeta Plus™ Depth Filters

3M™ Zeta Plus™ Depth Filters include a wide selection of single-use depth filters, featuring a full range of pore sizes and media types. The encapsulated filters deliver high contaminant capacity and throughput, eliminating a cleaning step and cleaning validation, and shorten setup in encapsulated systems.

[Learn more about the 3M™ Zeta Plus™ Depth Filters](#)



3M™ Emphaze™ AEX Hybrid Purifier

3M™ Emphaze™ AEX Hybrid Purifiers reduce harvest and clarification steps by replacing the second depth filter, accelerating operations and delivering high purity. Containing a Q-functional anion exchange media and a fine particle, bioburden reduction membrane, the 3M™ Emphaze™ consistently and efficiently reduces DNA, HCP, endotoxin levels and cell debris during manufacturing and downstream processing.

[Learn more about the 3M™ Emphaze™ AEX Hybrid Purifier](#)



3M™ Harvest RC Chromatographic Clarifier

3M™ Harvest RC Chromatographic Clarifier is the first one-step fiber chromatographic clarifier that accelerates the mAb harvest and clarification process by improving process efficiencies — like replacing the centrifuge and depth filters while delivering improved yield and purity.

[Learn more about the 3M™ Harvest RC Chromatographic Clarifier](#)



3M™ Polisher ST

3M™ Polisher ST is a fully encapsulated, scalable, single-use anion exchange (AEX) product with two complementary functionalized media for downstream processing. The AEX- functional membrane with novel ligand combined with a AEX-functional non-woven fiber allows it to perform in high salt conditions, delivers efficient removal of residual contaminants (HCP, DNA) and viral clearance even in the presence of turbidity.

[Learn more about the 3M™ Polisher ST](#)



3M™ LifeASSURE™ PDA Membrane Filters

3M™ LifeASSURE™ PDA Membrane Filters leverage single-use cartridges to help control bioburden and deliver process efficiencies through high capacity, fast flow rates and reliable filtration performance across gamma compatible and sterile formats.

[Learn more about the 3M™ LifeASSURE PDA Membrane Filters](#)



3M™ Encapsulated System Holders

3M™ Encapsulated System Holders are designed with operator ease and convenience in mind, accommodating a variety of 3M production capsules for biopharmaceutical processing lines. Multiple configurations of the holder designs allows for flexibility at various production scales, and provide small footprint during use and storage.

[Learn more about the 3M™ Encapsulated System Holders](#)



Take the next step with us.

Our service centers, quality management team and bioprocess application specialists are ready to work with you on your next challenge.

Visit https://www.3mindia.in/3M/en_IN/bioprocessing-in/
or contact your 3M representative to connect.



3M Purification Inc.
3M Separation and Purification Sciences Division
3M India Limited,
48-51, Electronic City, Phase-1,
Bangalore-560100

Phone 1800-425-3030
Web www.3mindia.in/3M/en_IN/bioprocessing-in/

Product Selection and Use: Many factors beyond 3M's control and uniquely within user's knowledge and control can affect the use and performance of a 3M product in a particular application. As a result, end-user is solely responsible for evaluating the product and determining whether it is appropriate and suitable for end-user's application, including completing a risk assessment that considers the product leachable characteristics and its impact on drug safety, conducting a workplace hazard assessment and reviewing all applicable regulations and standards. Failure to properly evaluate, select, and use a 3M product and appropriate safety products, or to meet all applicable safety regulations, may result in injury, sickness, death, and/or harm to property.

Warranty, Limited Remedy, and Disclaimer: Unless a different warranty is expressly identified on the applicable 3M product literature or packaging (in which case such express warranty governs), 3M warrants that each 3M product meets the applicable 3M product specification at the time 3M ships the product. 3M MAKES NO OTHER WARRANTIES OR CONDITIONS, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OR CONDITION OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR ARISING OUT OF A COURSE OF DEALING, CUSTOM, OR USAGE OF TRADE. If a 3M product does not conform to this warranty, then the sole and exclusive remedy is, at 3M's option, replacement of the 3M product or refund of the purchase price.

Limitation of Liability: Except for the limited remedy stated above, and except to the extent prohibited by law, 3M will not be liable for any loss or damage arising from or related to the 3M product, whether direct, indirect, special, incidental, or consequential (including, but not limited to, lost profits or business opportunity), regardless of the legal or equitable theory asserted, including, but not limited to, warranty, contract, negligence, or strict liability.

3M, Emphaze, LifeASSURE, and Zeta Plus are trademarks of 3M Company. All other trademarks are the property of their respective owners. © 3M 2022. All rights reserved. 70-2016-0360-5 REV 08/2022