

Self-Adherent Wrap

Non-Sterile

Product Description

3M™ Coban™ Self-Adherent Wrap is a cohesive elastic wrap constructed from a nonwoven material and elastic fibers. The cohesive properties allow the wrap to stick to itself but not to other materials or skin. The product is supplied non-sterile.

Indications for Use

Coban self-adherent wrap is intended for use as an elastic wrap to provide compression or support, or to secure dressings or devices. This product is not designed, sold or intended for use except as indicated.

Contraindications

None known.

Warnings

1. Coban self-adherent wrap will not slip or loosen after application. DO NOT wrap too tightly to avoid impairing circulation and causing serious injury. The risk of impaired circulation is higher for children.
2. Monitor area frequently (no less than every 24 hours) for signs of swelling, discoloration, pain, numbness, tingling, or other changes in sensation. If these symptoms occur, remove wrap immediately and contact your health care provider.
3. DO NOT place Coban self-adherent wrap in direct contact with an open wound because it may lead to infection.

Caution

This product contains natural rubber latex which may cause allergic reactions.

Precautions

1. When used as a component in a compression wrap system, it is important to ensure adequate arterial blood flow to avoid injury secondary to decreased circulation. Coban self-adherent wrap should be used as a component of a compression wrap system only under the supervision of a wound care specialist.
2. Do Not Reuse. Reuse may result in compromising product integrity and / or lead to device failure.

Directions for Use

1. Unwind approximately 12 inches of Coban self-adherent wrap from the roll and allow the wrap to relax.
2. Without stretching, apply one full wrap and overlap. Press overlapped area lightly to keep the end in place. See Figure 1.
3. Proceed with wrapping, applying tension to the wrap only to obtain desired compression. Always unwind more Coban self-adherent wrap from the roll before making a wrap so that the product is not applied with too much tension. See Figure 2.
4. When application is complete, cut off excess Coban self-adherent wrap. Secure the end of the application by gently pressing.

Storage/Shelf Life

This product should be stored at room temperature. Avoid excessive heat and humidity. For best performance, store Coban self-adherent wrap in its individual package until use. For shelf life, refer to the date printed on each primary carton.

How supplied

Supplied in boxes of individually packaged rolls. All rolls of Coban self-adherent wrap are 5 yards (4,5m) long stretched, except for 1584L which is 6.5 yards (6m) stretched.

If you have any questions or comments, please contact the 3M Health Care Customer Helpline at 1-800-228-3957.

Explanation of Symbols

 Caution: This product contains natural rubber latex which may cause allergic reactions.

 3M Health Care



St. Paul, MN 55144-1000

1-800-228-3957

www.3M.com/SkinWoundCare

© 2011, 3M. All rights reserved.

3M and Coban are trademarks of 3M.

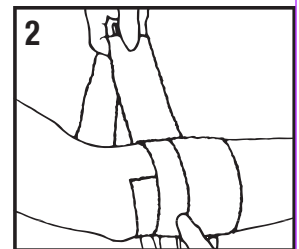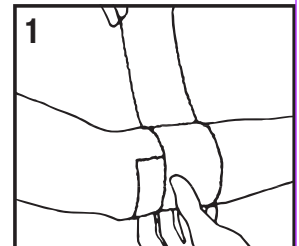

CE

EC REP

3M Health Care

D-41453 Neuss, Germany

Issue Date: 2011-03

34-8707-0789-9