

Product Description

3M™ Coban™ LF Latex Free Self-Adherent Wrap is a cohesive elastic wrap constructed from a nonwoven material and elastic fibers. The product does not contain natural rubber latex. The cohesive properties allow the wrap to stick to itself but not to other materials or skin. The product is supplied in individually packaged sterile rolls.

Indications for Use

Coban LF self-adherent wrap is intended for use as an elastic wrap to provide compression or support, or to secure dressings or devices. This product is not designed, sold or intended for use except as indicated.

Contraindications

None known

Warnings

1. Coban LF self-adherent wrap will not slip or loosen after application. DO NOT wrap too tightly to avoid impairing circulation and causing serious injury. The risk of impaired circulation is higher for children.
2. Monitor area frequently (no less than every 24 hours) for signs of swelling, discoloration, pain, numbness, tingling, or other changes in sensation. If these symptoms occur, remove wrap immediately and contact your Health Care provider.
3. DO NOT place Coban LF self-adherent wrap in direct contact with an open wound because it may lead to infection.

Precautions

1. When used as a component in a compression wrap system, it is important to ensure adequate arterial blood flow to avoid injury secondary to decreased circulation. Coban LF self-adherent wrap should be used as a component of a compression wrap system only under the supervision of a wound care specialist.
2. Do Not Reuse. Reuse may result in compromising product integrity and / or lead to device failure.

Directions for Use

1. Unwind approximately 12 inches of Coban LF self-adherent wrap from the roll and allow the wrap to relax.
2. Without stretching, apply one full wrap and overlap. Press overlapped area lightly to keep the end in place. See Figure 1.
3. Proceed with wrapping, applying tension to the wrap only to obtain desired compression. Always unwind more Coban LF self-adherent wrap from the roll before making a wrap so that the product is not applied with too much tension. See Figure 2.
4. When application is complete, tear off excess Coban self-adherent wrap. Secure the end of the application by gently pressing.

Storage/Shelf Life

This product should be stored at room temperature. Avoid excessive heat and humidity. For best performance, store Coban LF self-adherent wrap in its individual package until use. For shelf life, refer to the date printed on each package.

How supplied

Supplied in boxes of individually packaged rolls. All rolls of Coban LF self-adherent wrap are 5 yards (4,5m) long stretched. For sterile product, sterility is guaranteed unless individual package is damaged or opened.

If you have any questions or comments, please contact the 3M Health Care Customer Helpline at 1-800-228-3957.

Explanation of Symbols

 This product and package do not contain natural rubber latex


 Sterilized using irradiation

 Do not use if individual pouch is damaged or open

 Hand-tearable no scissors needed

 Caution, see instructions for use

 Do not reuse

 Use by date

 Batch code

 Manufacturer

 Stretched length

 Unstretched length




Made in U.S.A with globally sourced materials by

3M Health Care

2510 Conway Ave.
St. Paul, MN 55144
1-800-228-3957

© 3M, 2012. All rights reserved.
3M and Coban are trademarks of 3M.
www.3M.com/SkinWoundCare

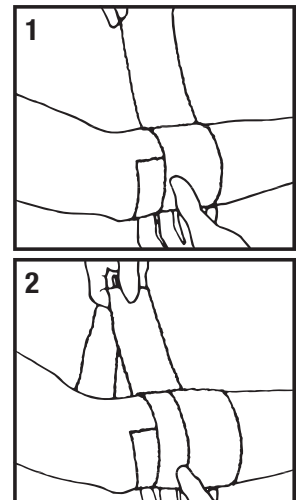

 **3M Deutschland GmbH**
Health Care Business
Carl-Schurz-Str. 1
41453 Neuss, Germany

Issue Date: 2012-01
34-8710-0048-4