


## 3M Infection Prevention Respiratory Protection

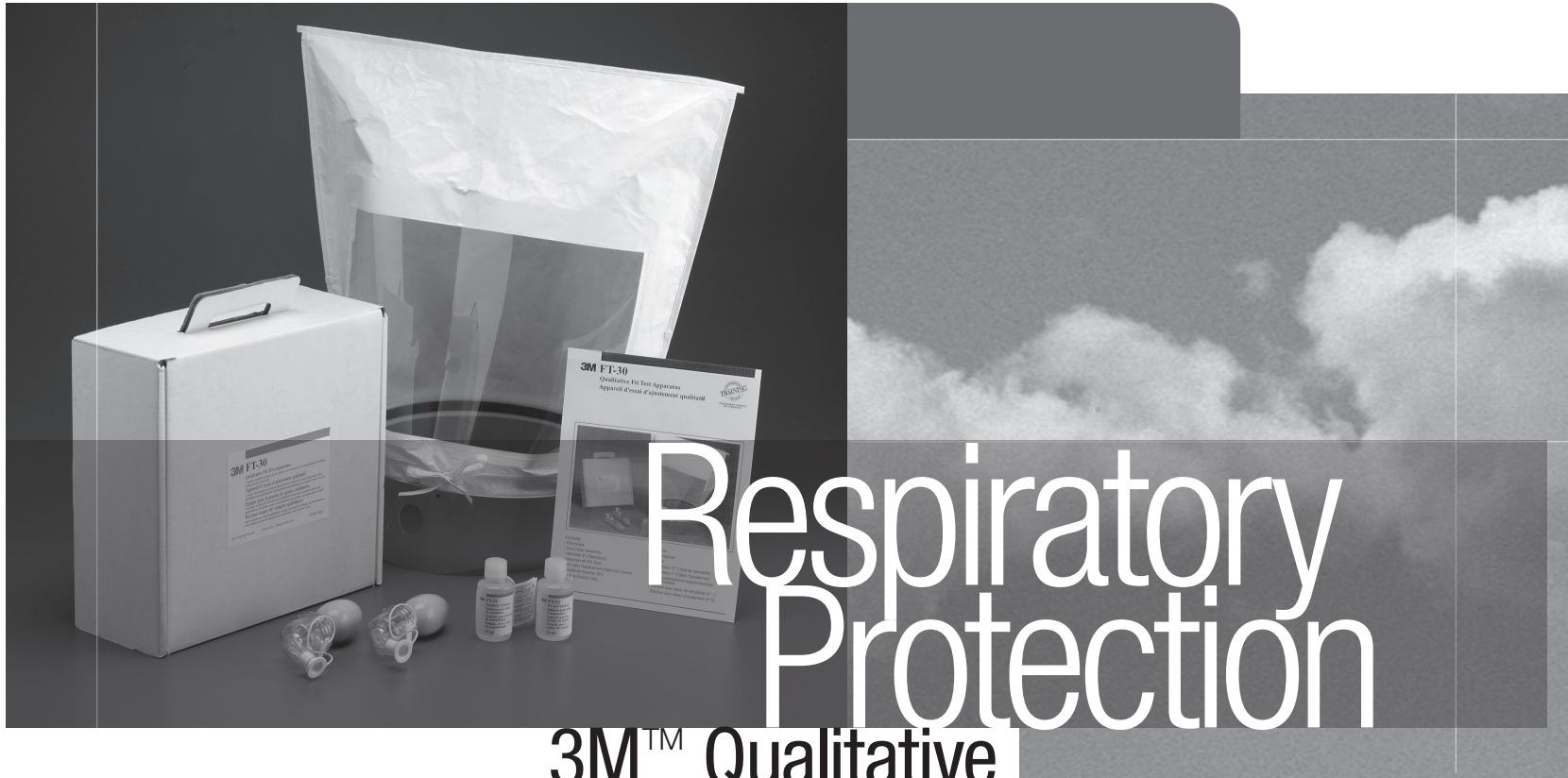

# Respiratory Protection

## 3M™ Qualitative Fit Test Apparatus

The 3M Qualitative Fit Test Apparatus offers a simple way to meet the Occupational Safety and Health Administration's (OSHA) qualitative fit testing requirements. It can be used to assess the face-to-respirator seal of any particulate respirator or gas/vapor respirator with a particulate prefilter.

Fit testing must be performed before a respirator is first issued to a healthcare worker and is required annually thereafter. (Note: a separate user seal check must be performed each time the respirator is worn.)

Fit testing is accomplished by spraying a test agent into the hood environment with the respirator in place. If the wearer does not detect the taste of the test agent, an acceptable fit has been demonstrated.












## 3M™ Qualitative Fit Test Apparatus



### Some benefits of 3M's Fit Test Apparatus:

- Fast, easy method for performing qualitative fit testing
- Two test agent options are available to determine employee's taste threshold to test solution – 3M™ Qualitative Fit Test Apparatus FT-10, Sweet and 3M™ Qualitative Fit Test Apparatus FT-30, Bitter
- Can be used with any particulate respirator, or any gas/vapor respirator with particulate prefilter
- Convenient, reliable and economical to use
- Provides consistent test atmosphere for test exercises
- Educational materials available from 3M to help you meet OSHA respirator program requirements
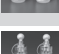
### Ordering Information

Cat. No.	Description	Size	Each/Box	Box/Case
 FT-10	3M™ Qualitative Fit Test Apparatus (sweet) Contains: FT-11, FT-12, FT-14, FT-15 (1 each); FT-13 (2 each)		1	1
 FT-30	3M™ Qualitative Fit Test Apparatus (bitter) Contains: FT-31, FT-32, FT-14, FT-15 (1 each); FT-13 (2 each)		1	1
 FT-11	Sensitivity Solution (sweet)	55 mL	1	6
 FT-12	Fit Test Solution (sweet)	55 mL	1	6
 FT-31	Sensitivity Solution (bitter)	55 mL	1	6
 FT-32	Fit Test Solution (bitter)	55 mL	1	6
 FT-13	Nebulizer		1	3
 FT-14	Test Hood (2/pack)		2	5
 FT-15	Collar		1	1

### Training, Fit Testing and Technical Service

Ask your 3M sales representative about our related training materials and services to help you meet OSHA training regulations for the use of respirators. Or call our technical service line at: 1-800-441-1922, opt. 3. Outside the U.S., contact your local 3M subsidiary.

For more information call the 3M Health Care Helpline at 1-800-228-3957 or visit our website at [www.3M.com/healthcare](http://www.3M.com/healthcare)



**3M Health Care**  
3M Center, Building 275-4W-02  
St. Paul, MN 55144-1000  
U.S.A.  
1 800 228-3957  
[www.3m.com/healthcare](http://www.3m.com/healthcare)

**3M Canada**  
Post Office Box 5757  
London, Ontario N6A 4T1  
Canada  
1 800 364-3577

3M is a trademark of 3M.  
Please recycle. Printed in U.S.A.  
© 3M 2008. All rights reserved.  
70-2008-6699-7