



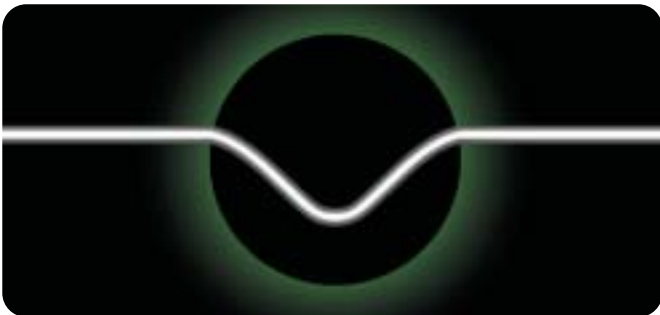

SMARTCLIP™  
SELF-LIGATING APPLIANCE SYSTEM

## DIMPLED ARCHWIRE

### A SIMPLE WAY TO AVOID “WALKING” ARCHWIRES

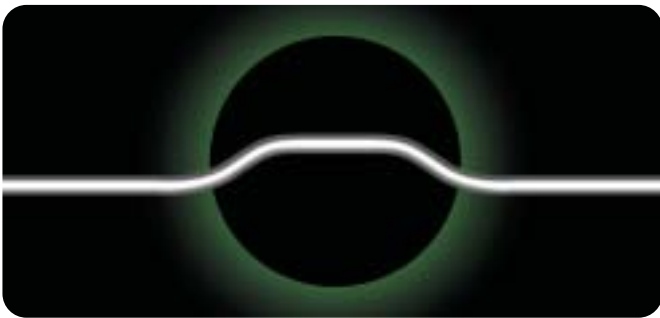

The shape and width of the dimpled archwires from 3M Unitek are arch-specific, to compensate for the differences in interbracket distances between the upper and lower teeth.



If you're among the growing number of orthodontists who have chosen the advantages of self-ligating appliances for patient treatment, or you are using minimal archwire ligation, then you (and just as importantly, your patients) may have experienced “walking” archwires.

Without ligation, archwires can move with relative ease along the brackets in the arch, and could reach the soft tissue at the back of the patient's mouth. This can lead to patient discomfort and possibly emergency appointments affecting you, your staff and the patient.

Dimpled archwires restrict archwire drift, thereby reducing the possibility of patient discomfort. Usable in all treatment scenarios, dimpled archwires reduce or eliminate the need to cinch the archwire behind the buccal tubes. The aesthetic dimple in the wires also gives a permanent visual mid-line reference for treatment.

#### DIMPLED ARCHWIRES FROM 3M UNITEK: AN EXCELLENT CHOICE

3M Unitek archwires are an exceptional choice for you and your patients. They're manufactured to the industry-acknowledged highest quality standards, and feature tight corner radii with an excellent surface finish. And, they're not just for self-ligating cases.

Choose from Nitinol Super-Elastic and Nitinol Classic wire, in five different sizes, round and rectangular, in the three OrthoForm™ Arch Forms. All are delivered in hygienic patient-specific unit-of-use packaging.

Take maximum advantage of the treatment opportunities of your self-ligating appliances, and reduce the opportunity for walking archwires. Choose SmartClip™ Appliance Dimpled Archwires for all your patients.



**SMARTCLIP™**  
SELF-LIGATING APPLIANCE SYSTEM

## DIMPLED ARCHWIRE


Wire sizes and types are easily identified with hygienic single-use, patient-specific packaging. Available in Unit of Use packaging only. 10 units per package.




upper




lower

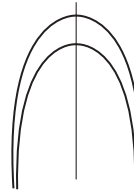

OrthoForm™ I  
Tapered

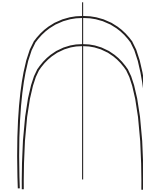

OrthoForm™ II  
Square

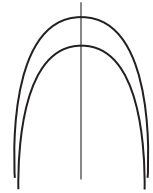

OrthoForm™ III  
Ovoid

Dimension	Alloy Type*	U/L	OrthoForm™ I Tapered	OrthoForm™ II Square	OrthoForm™ III Ovoid	
.014 ●	SE	U	9293-601	9296-805	9296-911	
	SE	L	9293-602	9296-806	9296-912	
	NC	U	9293-510	9296-405	9296-517	
	NC	L	9293-511	9296-406	9296-518	
	.016 ●	SE	U	9293-603	9296-807	9296-913
		SE	L	9293-604	9296-808	9296-914
NC		U	9293-512	9296-407	9296-519	
NC		L	9293-513	9296-408	9296-520	
.014 x .025 ■		SE	U	9293-813	9293-815	9293-817
		SE	L	9293-814	9293-816	9293-818
	NC	U	9293-819	9293-821	9293-823	
	NC	L	9293-820	9293-822	9293-824	
	.016 x .025 ■	SE	U	9293-849	9293-851	9293-853
		SE	L	9293-850	9293-852	9293-854
NC		U	9293-855	9293-857	9293-859	
NC		L	9293-856	9293-858	9293-860	
.017 x .025 ■		SE	U	9293-613	9297-807	9297-955
		SE	L	9293-614	9297-808	9297-956
	NC	U	9293-524	9297-707	9297-857	
	NC	L	9293-525	9297-708	9297-858	

\*NC (Nitinol Classic), SE (Nitinol Super-Elastic)

**3M** 3M Unitek  
Orthodontic Products

2724 South Peck Road  
Monrovia, CA 91016 USA

www.3MUnitek.com

In U.S. and Puerto Rico: 800-423-4588 • 626-574-4000

In Canada: 1-800-443-1661

Technical Helpline: 800-265-1943 • 626-574-4577

CE Hotline: 1-800-852-1990 x4649 • 626-574-4649

Outside these areas, contact your local representative.

©2006 3M Company.  
All rights reserved.  
016-945 0601