## Forsus<sup>™</sup> Fatigue Resistant Device EZ Module





## **Simply Click to Install**


The easy-to-install, "click-in-place" clip of the Forsus EZ Module is the latest innovation in the development of Class II Correctors. With the Forsus EZ Module, you now have an additional Forsus Corrector option of a fixed device that installs easily and resists breakage once engaged. The Forsus EZ Module provides a uniquely designed clip that "clicks" easily into a standard buccally-offset occlusal headgear tube by simply inserting it with any Weingart style plier. Mesial and distal stops on the clip assure a secure fit.


Its that easy. The Forsus EZ Module is an out-of-the-box solution that can improve your productivity and reduce uncomfortable patient chairtime with a simple "click." No "tweaks" or adjustments are usually required.

## **Patient Comfort and Treatment Efficiency**

The strength and security of the Forsus EZ Module is complemented by its design for patient comfort. An anti–rotational arm prevents the Forsus EZ Module from "rolling" and causing patient discomfort due to soft tissue irritation. The 3-piece telescoping spring fits perfectly with your current appliance system and the intraoral design is more attractive and comfortable for your patients.

Treatment efficiency is further enhanced by the tested fatigue resistance and consistent, long-lasting forces delivered by the Forsus EZ Module appliance. Tested fatigue resistance means less breakage, reducing the need to schedule spring-replacement appointments. Consistent forces yield more predictable outcomes in less time than you might require for other Class II Correctors, without the need for patient compliance. Other currently available solutions cannot compare with the Forsus EZ Module in terms of ease of installation, strength, reliability, and patient satisfaction.

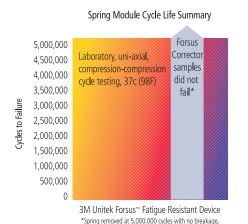

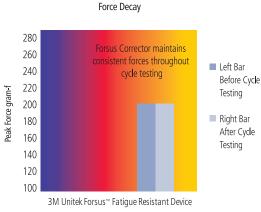

**3M** Unitek


Products that make your life easier.

Charts Source: 3M Unitek Laboratory

## Forsus<sup>™</sup> Fatigue Resistant Device EZ Module



Description	Unit		Part No.
Kits			
Forsus™ Fatigue Resistant Device EZ Module 20-Patient Kit	Kit		885-132
Forsus EZ Module 5-Patient Kit	Kit		885-131
ndividual Components			
Forsus EZ Module Reorder Pack (Engagement clip and spring)	pkg/5	L	885-138
Forsus EZ Module Reorder Pack (Engagement clip and spring)	pkg/5	R	885-139
Direct Push Rod Short (25mm)	pkg/5	L	885-111
Direct Push Rod Short (25mm)	pkg/5	R	885-112
Direct Push Rod Medium (29mm)	pkg/5	L	885-113
Direct Push Rod Medium (29mm)	pkg/5	R	885-114
Direct Push Rod Large (32mm)	pkg/5	L	885-115
Direct Push Rod Large (32mm)	pkg/5	R	885-116
Direct Push Rod X-Large (35mm)	pkg/5	L	885-117
Direct Push Rod X-Large (35mm)	pkg/5	R	885-118
Direct Push Rod XX-Large (38mm)	pkg/5	L	885-107
Direct Push Rod XX-Large (38mm)	pkg/5	R	885-108
Jniversal Split Crimp	pkg/50		885-203
Measurement Guide	pkg/5		807-014

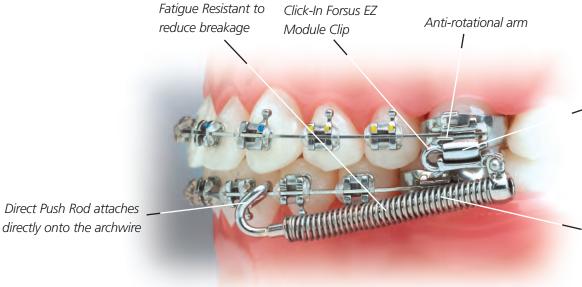

Buccally-offset occlusal headgear tubes fit securely into the clip channel

A 3-piece telescoping coaxial spring assembly designed to allow for enhanced patient jaw opening

