

Victory Series™ Funneled 2nd Molar Buccal Tubes


For Orthodontists concerned about root torque on their patient's second molars, Victory Series™ Funneled 2nd Molar Tubes are an excellent choice. These tubes offer a range of torques, specifically with 2nd molar torque control in mind. Victory Series 2nd molar tubes address treatment issues that can arise due to lingual cusp interference during the finishing stages of treatment. The funneled design also helps ease archwire insertion.

Works Great With 3M™ Self-Ligating Appliance Systems

These tubes are ideal to use with SmartClip™ and Clarity™ SL Self-Ligating Appliance Systems, as archwire engagement is enhanced with the wide opening entryway.

Victory Series 2nd Molar Tube Features


- Higher torque values acknowledge that upper 2nd molars have different tooth morphology than upper 1st molars. The higher torque helps to better align the upper 2nd molar lingual cusp with the rest of the arch.
- The kind of design quality expected of 3M and its legacy of quality buccal tubes: crafted from metal injection molding technology, which yields
 - Consistent quality and dimension accuracy
 - Contoured design for patient comfort
 - Strength in a low-profile design
- Fully compatible with Victory Series, SmartClip SL3 and Clarity SL appliance systems
- Funneled mesial entry to assist with wire insertion
- Contoured base design for bond strength and facilitated fit
- Adjustable low profile hook for versatile function and patient comfort
- Side gripping notches for ease of handling and placement
- Distal notch for ease of ligation
- Alignment seating notch to assist with more accurate placement
- Available with APC™ II and APC™ PLUS Adhesive Coating

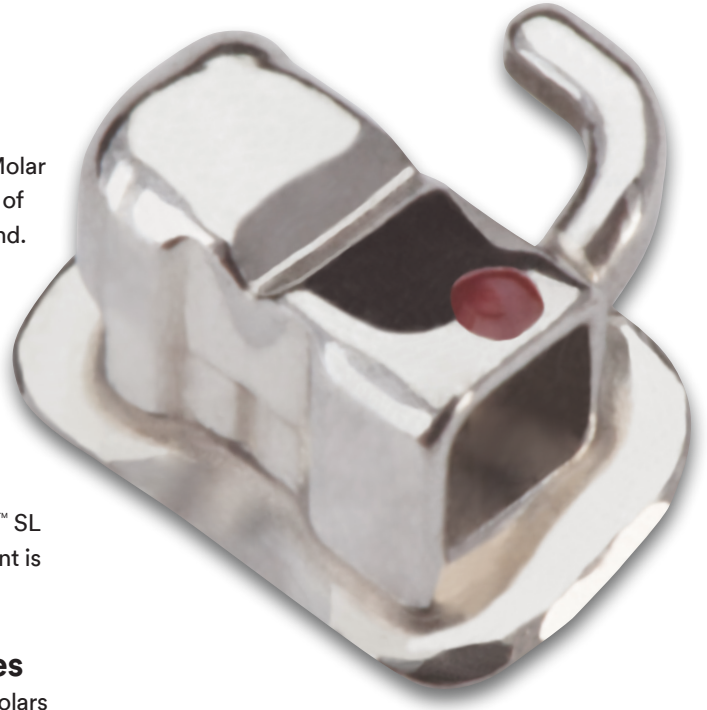

The Victory Series 2nd molar tubes exemplify the qualities of efficiency and functionality that are at the core of the MBT™ Versatile+ Appliance System and Sondhi System treatment philosophies.

Victory Series™

Funneled 2nd Molar Buccal Tubes

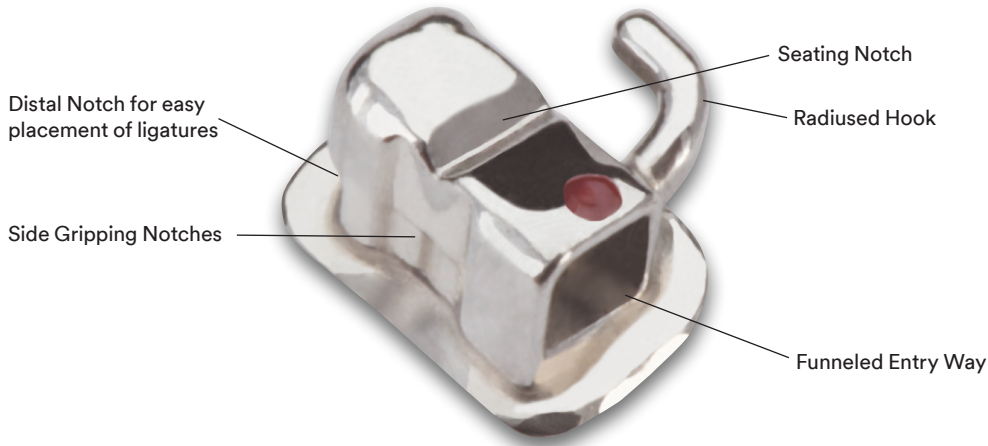

Upper



Lower

Part Numbers and Specifications

Victory Series™ Non-Convertible 2nd Molar Buccal Tubes

Maxillary										
Slot Torque	Distal Offset	Length	Slot	L/R	APC™ PLUS	APC™ II	Part Number	Base	System	
-19°	8°	3.9 mm	.018 in	UL	5	3	067-873	1	MBT™ / Sondhi	
-19°	8°	3.9 mm	.018 in	UR	5	3	067-872	2	MBT / Sondhi	
-19°	8°	3.9 mm	.022 in	UL	5	3	067-973	1	MBT / Sondhi	
-19°	8°	3.9 mm	.022 in	UR	5	3	067-972	2	MBT / Sondhi	
Mandibular										
Slot Torque	Distal Offset	Length	Slot	L/R	APC™ PLUS	APC™ II	Part Number	Base	System	
-10°	0°	3.6 mm	.018 in	LL	5	3	067-536	3	MBT	
-10°	0°	3.6 mm	.018 in	LR	5	3	067-537	4	MBT	
-10°	0°	3.6 mm	.022 in	LL	5	3	067-636	3	MBT	
-10°	0°	3.6 mm	.022 in	LR	5	3	067-637	4	MBT	
-25°	6°	3.6 mm	.018 in	LL	5	3	067-538	3	Sondhi	
-25°	6°	3.6 mm	.018 in	LR	5	3	067-539	4	Sondhi	
-25°	6°	3.6 mm	.022 in	LL	5	3	067-638	3	Sondhi	
-25°	6°	3.6 mm	.022 in	LR	5	3	067-639	4	Sondhi	

APC™ PLUS

APC™ II

Base

Use these prefixes to order APC™ PLUS or APC™ II Adhesive Coated Appliance System.

APC™ Adhesive coating only available with base.



3M Oral Care
2510 Conway Avenue
St. Paul, MN 55144-1000 USA
Phone 1-800-423-4588
Web 3M.com/ortho

3M Canada
Health Care Division
300 Tartan Dr.
London, ON N5V 4M9 Canada
Phone 1-800-443-1661

3M, APC, Clarity, MBT, SmartClip, and Victory Series are trademarks of 3M. Used under license in Canada. © 3M 2016. All rights reserved. 70-2021-3942-7 1601