

Tegaderm™ Dressing Selection for the Moisture Management Continuum



LOW

Managing **Low to Moderate** Levels of Exudate

Tegaderm™ Absorbent Clear Acrylic

- Ideal for low to moderately exuding wounds such as skin tears, abrasions, pressure ulcers, surgical wounds and donor sites
- Clear dressing allowing wound bed to be monitored without dressing removal
- Designed for longer wear time - needs to be changed if dressing is leaking, lifting off or there is wound fluid under the adhesive border



Wound Exudate Level

MOD

Managing **Moderate to High** Levels of Exudate



Tegaderm™ High Performance Foam Non-Adhesive

- Ideal for moderate to highly exuding wounds such as pressure ulcers, venous and arterial ulcers
- Non-adhesive foam can be secured with tape, compression wraps, gauze, self-adherent wrap or other methods
- Ideal choice for use with compression or when all adhesives should be avoided



Tegaderm™ Silicone Foam

- Ideal for moderate to highly exuding wounds such as pressure ulcers, venous and arterial ulcers
- Stays in place with a gentle 3M™ Silicone Adhesive and provides atraumatic removal
- Ideal choice for patients with at-risk or fragile skin or if more frequent dressing changes are needed



Tegaderm™ High Performance Foam Adhesive

- Ideal for moderate to highly exuding wounds such as pressure ulcers, venous and arterial ulcers
- Ideal for most cases, especially where increased securement or longer wear time is desired
- Secures best to challenging areas such as heels and elbows, and to both moist and dry skin
- Designed for wear time of up to seven days

HIGH

To avoid all adhesives on the skin

Fragile skin
Frequent dressing changes expected

Ideal for most cases
Most reliable wear time

Tegaderm™ Dressing Selection for the Moisture Management Continuum




Tegaderm™ Foam Border Dressings

The science of securing to the skin: two adhesive choices offer the flexibility you need


The Tegaderm™ Foam Family of Dressings offers clinicians two adhesive systems to meet your patient needs. Tegaderm™ High Performance Foam Adhesive Dressing featuring an acrylate adhesive, is ideal for most cases as a general use foam dressing, securing especially well to difficult body contours (heels, elbows, etc.) as well as to both dry and moist skin. Tegaderm™ Silicone Foam Border Dressing is designed with a gentle 3M™ Silicone Adhesive and is the ideal choice for patients with fragile and at-risk skin.

Acrylate adhesive

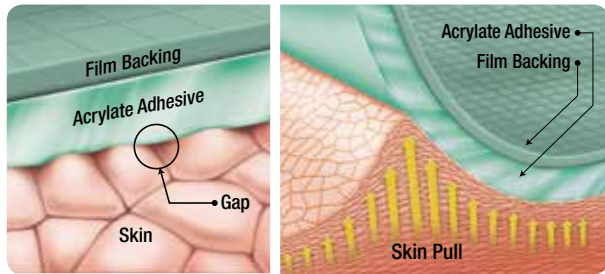

Acrylate adhesive initially adheres to the skin cells closest to the top, leaving some gaps in adherence. Over time the adhesive fills the gaps and strengthens as it forms a tighter bond with the skin, making it ideal for situations where increased securement or longer wear are needed.

Silicone adhesive

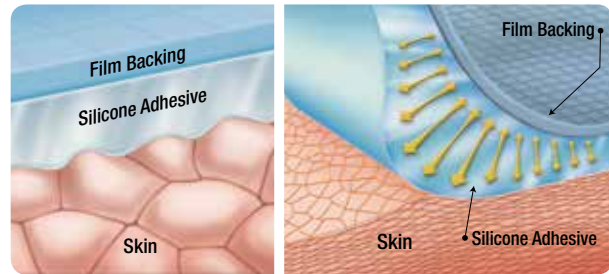

Silicone adhesive conforms to the uneven surfaces of the skin immediately and remains constant throughout wear time. Unlike acrylate adhesives, silicone adhesives have lower surface tension, making them the preferred choice for those patients with at-risk or fragile skin or when more frequent dressing changes are required.

Product Ordering Guide – Moisture Management Dressings

3M Cat. ID	Product Description	3M Cat. ID	Product Description
Tegaderm™ Absorbent Clear Acrylic Dressing		Tegaderm™ High Performance Foam Non-Adhesive Dressing	
90800	Oval, 3.8 cm x 5.7 cm (pad), 7.6 cm x 9.5 cm (overall)	90600	5 cm x 5 cm
90801	Oval, 6.0 cm x 7.6 cm (pad), 11.1 cm x 12.7 cm (overall)	90601	10 cm x 10 cm
90802	Square, 9.8 cm x 10.1 cm (pad), 14.9 cm x 15.2 cm (overall)	90602	10 cm x 20 cm
90803	Oval, 8.5 cm x 10.7 cm (pad), 14.2 cm x 15.8 cm (overall)	90603	20 cm x 20 cm
90805	Square, 14.9 cm x 15.2 cm (pad), 20 cm x 20.3 cm (overall)	90604	Fenestrated, 8.9 cm x 8.9 cm
90807	Sacral, 11.75 cm x 14.3 cm (pad), 16.8 cm x 19 cm (overall)	90605	Roll, 10 cm x 60 cm
Tegaderm™ High Performance Foam Adhesive Dressing		Tegaderm™ Silicone Foam Border Dressing	
90610	Square, 5 cm x 5 cm (pad), 8.9 cm x 8.9 cm (overall)	90640	Square, 3.7 cm x 3.7 cm (pad), 7.5 cm x 7.5 cm (overall)
90611	Oval, 6 cm x 7.6 cm (pad), 10 cm x 11 cm (overall)	90641	Square, 5.9 cm x 5.9 cm (pad), 10 cm x 10 cm (overall)
90612	Square, 10 cm x 10 cm (pad), 14.3 cm x 14.3 cm (overall)	90642	Square, 10.3 cm x 10.3 cm (pad), 15 cm x 15 cm (overall)
90613	Oval, 10 cm x 11 cm (pad), 14.3 cm x 15.6 cm (overall)		
90614	Oval, 3.1 cm x 3.8 cm (pad), 6.9 cm x 7.6 cm (overall)		
90616	Oval, 14 cm x 17.1 cm (pad), 19 cm x 22.2 cm (overall)		
90619	Heel/Elbow, 7.6 cm x 7.6 cm (pad), 14 cm x 14 cm (overall)		



Available in Canada from:
Critical & Chronic Care Solutions Division
3M Canada Company
P.O. Box 5757
London, Ontario N6A 4T1
Canada
1 800-364-3577
www.3M.com

3M Health Care
2510 Conway Avenue
St. Paul, MN 55144-1000
U.S.A.
1 800-228-3957
www.3M.com/healthcare

3M and Tegaderm are trademarks of 3M.
Used under license in Canada.
© 2014, 3M. All rights reserved.
Please recycle. Printed in Canada.
1405-01052E