

3M post and core solution. Complexity simplified.

For predictable results – in less time: The endodontic procedure is very complex. So the simpler and more reliable the endodontic restorative procedure is, the better – for both the dentist and the patient. With fewer clinical steps and fewer products, the 3M post and core solution lets you work faster, easier and more efficiently.

And, it helps you deliver predictable, robust outcomes. All you need is the new 3M™ RelyX™ Fiber Post 3D Glass Fiber Post plus three other easy-to-use 3M products: 3M™ RelyX™ Unicem 2 Automix Cement, 3M™ Scotchbond™ Universal Adhesive, and 3M™ Filtek™ One Bulk Fill Restorative. These four products have one common objective: Simplify the complex!



3M™ Filtek™ One
Bulk Fill Restorative




3M™ Scotchbond™
Universal
Adhesive



3M™ RelyX™ Fiber Post 3D
Glass Fiber Post



3M™ RelyX™ Unicem 2 Automix
Self-Adhesive Universal
Resin Cement

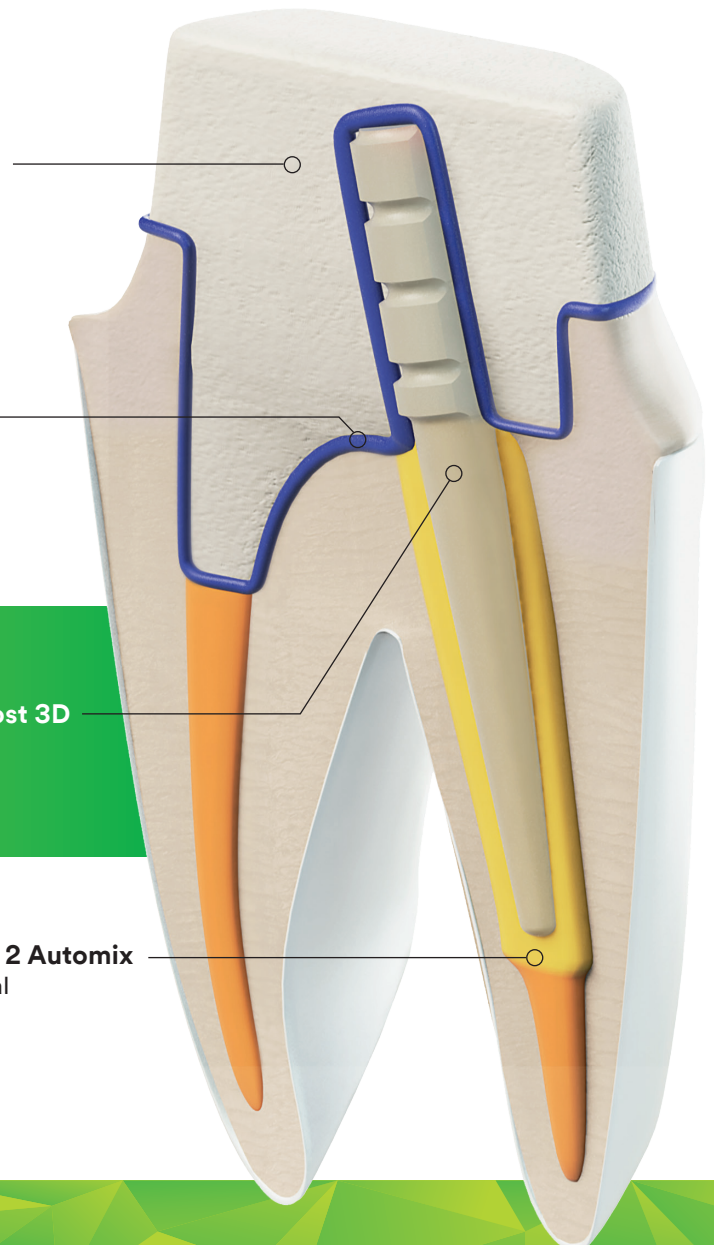

3M™ RelyX™ Unicem 2 Automix Self-Adhesive Universal Resin Cement offers excellent bond strength plus saves you precious time, as there is no need to etch, prime or bond. The automix syringe, together with the endo tip, lets you apply the ready-mixed cement directly in the root canal from bottom to top. In one application – without voids.

3M™ RelyX™ Fiber Post 3D Glass Fiber Post's coronal 3D macro retention offers high mechanical retention to the core build-up material. Due to its microporous surface it also provides a strong bond in the root canal without a silanating step. Post position control is greatly improved as a result of its enhanced radiopacity.

3M™ Scotchbond™ Universal Adhesive provides a reliable bond of the core build-up material to the post and tooth structure – with no need for an additional primer. One application – done!

3M™ Filtek™ One Bulk Fill Restorative has excellent material and handling characteristics that make Filtek One Bulk Fill Restorative a great choice of core build-up material. With one-step increments up to 5mm, core build-ups can be done faster and easier.



Ordering Information

Item No.	Product Information
3M™ RelyX™ Fiber Post 3D Glass Fiber Post Kits	
56957	Introductory Kit: 20 posts sizes 0, 1, 2, 3 (5 posts of each size), 5 drills sizes 0, 1, 2, 3, Universal (1 drill of each size)
56958	Trial Kit with RelyX™ Unicem 2 Automix Cement: 10 posts sizes 1, 2 (5 posts of each size), 3 drills sizes 1, 2, Universal (1 drill of each size), 1 – 3M™ RelyX™ Unicem 2 Automix Syringe A2 – 8.5g, 10 mixing tips, 10 endo tips
56960	Post and Core Build-up Kit with RelyX™ Unicem 2 Automix Cement: 1 Trial Kit 3M™ RelyX™ Fiber Post 3D/ 3M™ RelyX™ Unicem 2 Automix Cement, 20 capsules 3M™ Filtek™ One Bulk Fill Restorative A3 – 0.2g each, 1 vial 3M™ Scotchbond™ Universal Adhesive – 5ml
3M™ RelyX™ Fiber Post 3D Glass Fiber Posts and Drills	
56948	10 posts, size 0 (1.1 mm dia), white
56949	10 posts, size 1 (1.3 mm dia), yellow
56950	10 posts, size 2 (1.6 mm dia), red
56951	10 posts, size 3 (1.9 mm dia), blue
56873	1 drill, size 0 (1.1 mm dia), white
56864	1 drill, size 1 (1.3 mm dia), yellow
56865	1 drill, size 2 (1.6 mm dia), red
56866	1 drill, size 3 (1.9 mm dia), blue
56867	1 Universal drill
3M™ RelyX™ Unicem 2 Automix Self-Adhesive Resin Cement	
	Refills – each containing 1 syringe – 8.5g, 10 mixing tips (regular), 5 mixing tips (wide), 5 endo tips
56846	Shade A2 Universal
56847	Shade A3 Opaque
56848	Shade Translucent
	Value Packs – each containing 3 syringes – 8.5g, 30 mixing tips (regular), 20 mixing tips (wide), 15 endo tips, 5 intra-oral tips
56849	Shade A2 Universal
56857	Shade A3 Opaque
56858	Shade Translucent
Accessories – Cement Mixing Tips	
56851	30 mixing tips (regular)
56919	15 mixing tips (wide)
56918	15 mixing tips (wide), 15 intra-oral tips

Item No.	Product Information
3M™ Scotchbond™ Universal Adhesive	
41254	Introductory Kit: 1 vial 3M™ Scotchbond™ Universal Adhesive – 5 ml, 50 disposable mixing wells, 50 disposable microbrushes, 1 etchant syringe – 3ml, 25 dispensing tips for etchant syringe
41258	Refill: 1 vial 3M™ Scotchbond™ Universal Adhesive – 5 ml
3M™ Filtek™ One Bulk Fill Restorative	
Refill	Item No. Item No.
Shade	Syringe incl. 1 syringe – 4g Capsule incl. 20 capsules – 0.2g each
A1	4866A1 4867A1
A2	4866A2 4867A2
A3	4866A3 4867A3
B1	4866B1 4867B1
C2	4866C2 4867C2
3M™ Filtek™ Bulk Fill Flowable Restorative	
Refill	Item No. Item No.
Shade	Syringe incl. 2 syringes – 2g each, 20 dispensing tips Capsule incl. 15 capsules – 0.2g each
A1	4862A1 4861A1
A2	4862A2 4861A2
A3	4862A3 4861A3
Universal	4862U 4861U
3740T	20 dispensing tips 19-gauge



3M Oral Care
2510 Conway Avenue
St. Paul, MN 55144-1000
USA
1-800-634-2249

3M Canada
Post Office Box 5757
London, Ontario N6A 4T1
Canada
1-888-363-3685

3M.com/RelyX3D

3M, ESPE, Filtek, RelyX and Scotchbond are trademarks of 3M Company or 3M Deutschland GmbH. Used under license in Canada. All other trademarks are owned by other companies. Please recycle. Printed in U.S.A. © 3M 2017. All rights reserved.
70-2013-0726-4 (Rev. 9/2017)