

3M

Filtek™

**One
Bulk Fill Restorative**



**You choose it
for speed.
Now it's more
esthetic, too.**



You'll like that it's simple.
They'll be happy that it
looks good, too.

3M™ Filtek™ One Bulk Fill Restorative



It's true. Now bulk fill placement can be more esthetic.

Posterior restorations are challenging—which is why bulk fill composites that simplify and speed up the procedure have become so popular.


However, until now, the faster, easier placement of bulk fill has been at the sacrifice of esthetics. So, as a global leader in esthetic and restorative technologies, we went back to our experts for a solution you're sure to love:

3M™ Filtek™ One Bulk Fill Restorative.

Designed for the posterior, its unique optical properties and improved opacity give you the simplicity of one-step placement—up to 5 mm—*without compromising esthetic results.*

Now you and your patients have something to smile about.

- Fast and easy one-step placement
- Increased opacity for improved esthetics
- Stress relief to enable up to 5mm depth of cure
- Excellent adaptation
- No need for layering or expensive dispensing devices
- Excellent handling and sculptability
- 3M's true nanotechnology provides superior wear resistance and excellent polish retention
- High radiopacity
- Unique tip design for easier access to deep cavities
- Available in 5 shades

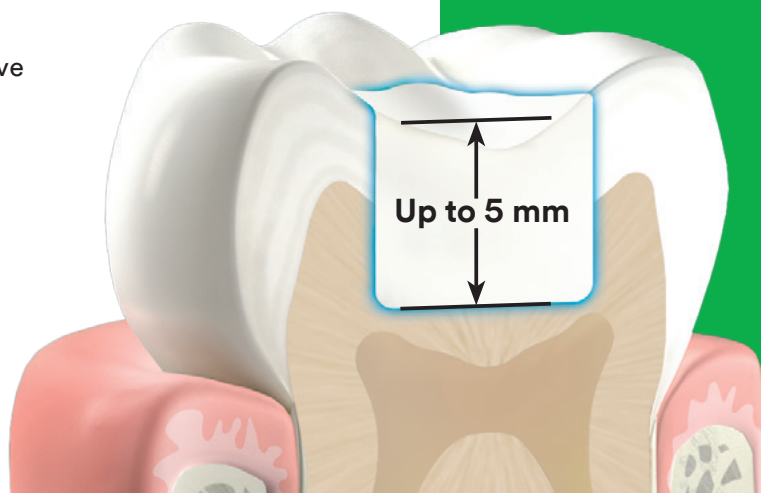

How 3M Science enables the simplicity of one-step bulk placement with improved esthetics.

The technologies behind 3M™ Filtek™ One Bulk Fill Restorative

NEW! Smart Contrast Ratio Management

Increases opacity without sacrificing depth of cure


Further development of 3M nanotechnology allows us to manage opacity and depth of cure in a smarter way. With our unique Smart Contrast Ratio Management technology, we've increased opacity without reducing depth of cure.

In its cured state, 3M™ Filtek™ One Bulk Fill Restorative has a higher opacity than 3M™ Filtek™ Bulk Fill Posterior Restorative, resulting in improved esthetics.

Did you know?

Contrast Ratio is the measurement of the translucency or opacity of the material.


Contrast Ratio

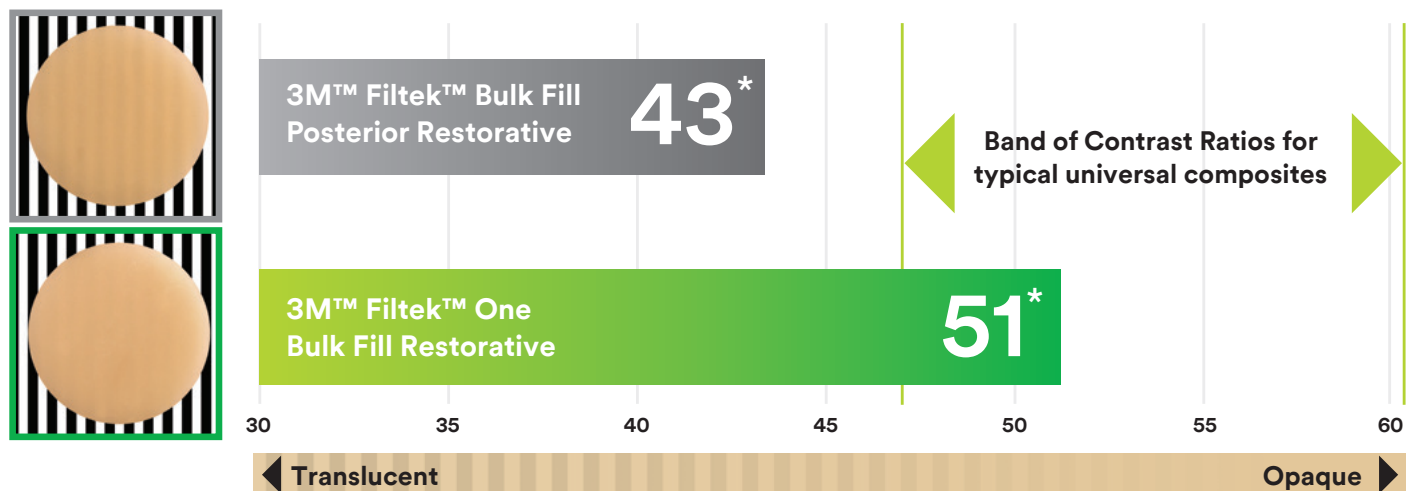


*The contrast ratio is the average of all shades.

Innovative AUDMA and AFM monomers

Reduces shrinkage and relieves stress to allow bulk fill placement up to 5 mm.

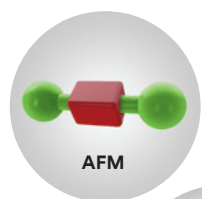

Addition-fragmentation monomer (AFM) relieves stress during polymerization.

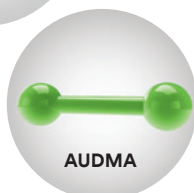


Aromatic urethane dimethacrylate (AUDMA) helps reduce the amount of shrinkage and stress that occurs during polymerization.


3M's TRUE nanofiller technology

Provides superior wear resistance and excellent polish retention.

5K magnification of nanocluster

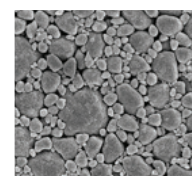

50K magnification of nanocluster

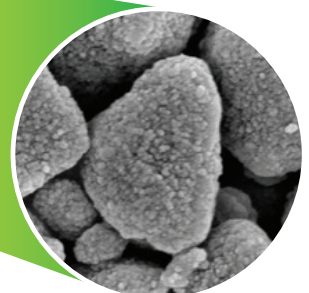

Now that you've repaired the tooth, help them protect the restoration at home.

3M™ Clinpro™ Tooth Crème 0.21% Sodium Fluoride Anti-Cavity Toothpaste



Contains 0.21% NaF (950 ppm)

Non-prescription-strength formula 3M™ Clinpro™ Tooth Crème 0.21% Sodium Fluoride Anti-Cavity Toothpaste is the answer. In addition to a lower level of fluoride, it contains functionalized Tri-Calcium Phosphate (FTCP), available exclusively from 3M—so it releases calcium and phosphate at the tooth surface, which helps to maintain healthy teeth.

- ✓ Protects against cavities and helps reverse white spots
- ✓ Proven through 5 years of clinical use
- ✓ Cleans and whitens teeth with low abrasion
- ✓ Replaces the patient's current toothpaste in the toothbrushing step, making compliance easy

Ordering Information



3M™ Filtek™ One Bulk Fill Restorative

Refill	Item #	Item #
Shade	Syringe	Capsule
A1	4866A1	4867A1
A2	4866A2	4867A2
A3	4866A3	4867A3
B1	4866B1	4867B1
C2	4866C2	4867C2

SYRINGE SHADE REFILLS include: 1 – 4g syringe; Technique Guide; Instructions for Use.

CAPSULE SHADE REFILLS include: 20 – 0.2g capsules; Technique Guide; Instructions for Use.



3M™ Clinpro™ Tooth Crème 0.21% Sodium Fluoride Anti-Cavity Toothpaste

Item #	Product Information
12218	Vanilla Mint Flavor – 4 oz tube
12218SM	Spearmint Flavor – 4 oz tube

www.3M.com



3M Oral Care
2510 Conway Avenue
St. Paul, MN 55144-1000 USA
Phone 1-800-634-2249
Web 3M.com/dental

3M Deutschland GmbH
Location Seefeld
3M ESPE • ESPE Platz
82229 Seefeld • Germany
Info3MESPE@mmm.com
www.3MESPE.com

3M, ESPE, Clinpro and Filtek are trademarks of 3M or 3M Deutschland GmbH. Used under license in Canada.
© 3M 2016. All rights reserved.