

3M post and core solution. Complexity simplified.

For predictable results – in less time: The endodontic procedure is very complex. So the simpler and more reliable the endodontic restorative procedure is, the better – for both the dentist and the patient. With fewer clinical steps and fewer products, the 3M post and core solution lets you work faster, easier and more efficiently.

And, it helps you deliver predictable, robust outcomes. All you need is the new RelyX™ Fiber Post 3D Glass Fiber Post plus three other easy-to-use 3M products: RelyX™ Unicem 2 Automix Cement, Scotchbond™ Universal Adhesive, and Filtek™ Bulk Fill Posterior Restorative. These four products have one common objective: Simplify the complex!



Filtek™ Bulk Fill
Posterior Restorative




Scotchbond™ Universal
Adhesive



RelyX™ Fiber Post 3D
Glass Fiber Post



RelyX™ Unicem 2 Automix
Self-Adhesive Resin Cement

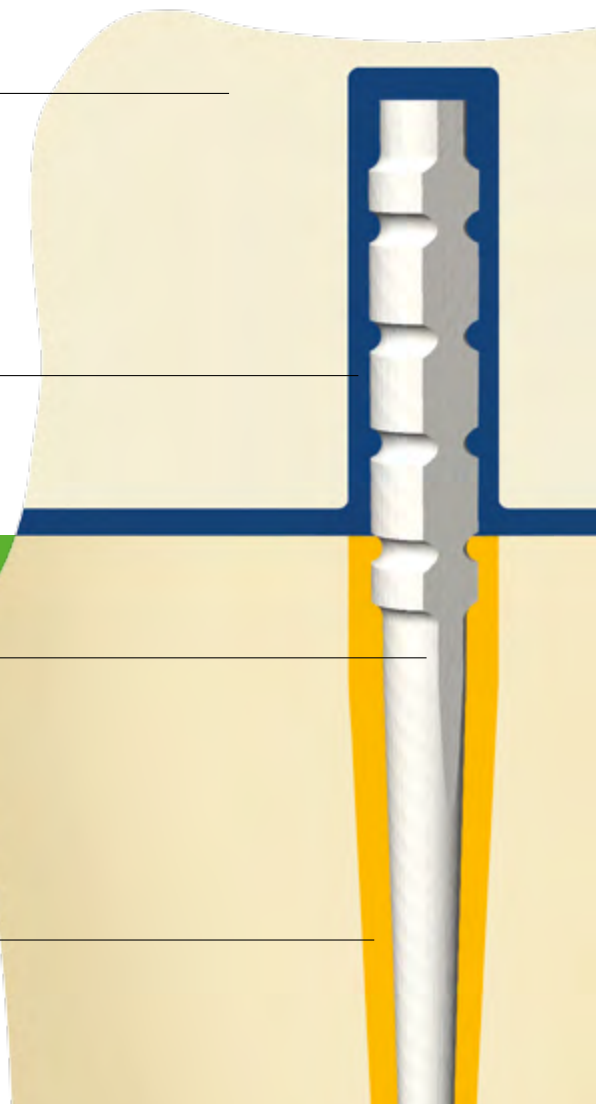

RelyX™ Unicem 2 Automix Self-Adhesive Resin Cement offers excellent bond strength plus saves you precious time, as there is no need to etch, prime or bond. The automix syringe, together with the endo tip, lets you apply the ready-mixed cement directly in the root canal from bottom to top. In one application – without voids.

RelyX™ Fiber Post 3D Glass Fiber Post's coronal 3D macro retentions offer high mechanical retention to the core build-up material. Due to its microporous surface it also provides a strong bond in the root canal without a silanating step. Post position control is greatly improved as a result of its enhanced radiopacity.

Scotchbond™ Universal Adhesive provides a reliable bond of the core build-up material to the post and tooth structure – with no need for an additional primer. One application – done!

Filtek™ Bulk Fill Posterior Restorative is ideally suited as core build-up material, especially for big lesions in combination with a post: The capsule design ensures a void-free application. The material is sculptable but also easy to adapt to both tooth structure and post. And, for easier preparation, the cured material cuts like dentine. Core build-ups can be done quickly and easily.



Ordering Information

Item No.	Product Information
RelyX™ Fiber Post 3D Glass Fiber Post	
56958	Trial Kit with RelyX™ Unicem 2 Automix Cement: 10 posts sizes 1, 2 (5 posts of each size), 3 drills sizes 1, 2, Universal (1 drill of each size), 1 RelyX™ Unicem 2 Automix Syringe A2 – 8.5g, 10 mixing tips, 10 endo tips
Refills – posts and drills	
56948	10 posts, size 0 (1.1 mm dia), white
56949	10 posts, size 1 (1.3 mm dia), yellow
56950	10 posts, size 2 (1.6 mm dia), red
56951	10 posts, size 3 (1.9 mm dia), blue
56873	1 drill, size 0 (1.1 mm dia), white
56864	1 drill, size 1 (1.3 mm dia), yellow
56865	1 drill, size 2 (1.6 mm dia), red
56866	1 drill, size 3 (1.9 mm dia), blue
56867	1 Universal drill
RelyX™ Unicem 2 Automix Self-Adhesive Resin Cement Refills – each containing 1 syringe – 8.5g, 15 mixing tips (regular), 5 mixing tips (wide), 5 endo tips	
56846	Shade A2 Universal
56847	Shade A3 Opaque
56848	Shade Translucent
Mixing tips	
56851	Mixing tips - regular (30)
56918	Mixing tips - wide (15); intra-oral tips (15)
56919	Mixing tips - wide (15); endo tips (15)

Item No.	Product Information	
Scotchbond™ Universal Adhesive		
41254	Introductory Kit: 1 vial Scotchbond™ Universal Adhesive – 5 ml, 50 disposable mixing wells, 50 disposable microbrushes, 1 etchant syringe – 3 ml, 25 dispensing tips for etchant syringe	
41258	Refill: 1 vial Scotchbond™ Universal Adhesive – 5 ml	
Filtek™ Bulk Fill Posterior Restorative		
Refill	Item No.	Item No.
Shade	Syringe incl. 1 syringe – 4g	Capsule incl. 20 capsules – 0.2g each
A1	4863A1	4864A1
A2	4863A2	4864A2
A3	4863A3	4864A3
B1	4863B1	4864B1
C2	4863C2	4864C2



3M Oral Care

3M Australia Pty Limited
 Building A, 1 Rivett Rd
 North Ryde, NSW 2113
 Ph: 1300 363 454
 3M.com.au/dental

3M New Zealand Limited
 94 Apollo Drive
 Rosedale, Auckland 0632
 Ph: 0800 80 81 82
 3M.co.nz/dental

Always follow the Instructions For Use (IFU) and refer to IFU for full indications, precautions and warnings. Claims supported by 3M Data on file. Information correct at time of printing. 3M, “3M Science. Applied to Life.”, Filtek, RelyX and Scotchbond are trademarks of 3M. Please recycle. © 3M 2017. All rights reserved.