

3M Science.
Applied to Life.™

3M™ Impregum™ Super Quick Polyether Impression Material

Capture every
detail in 2 minutes?
Yes, it can.



Quicker set.

With just a two minute setting time the new Impregum Super Quick Material allows for faster impressions in everyday cases like single-unit crowns, implants, inlays/onlays and smaller bridges up to 3 units.



100% Polyether. Trusted precision.

Like all polyethers Impregum Super Quick Materials demonstrate excellent performance in moist situations and exhibit great flow behaviour, leading to higher levels of accuracy and detailed, high-precision restorations.

The nature of polyether allows you to get the full 45 seconds of working time, avoiding inaccuracies associated with premature setting times.



Make your patients feel comfortable.

An improved fresh, minty taste together with a short time in the mouth significantly increases your patients' comfort.

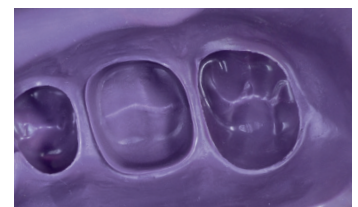
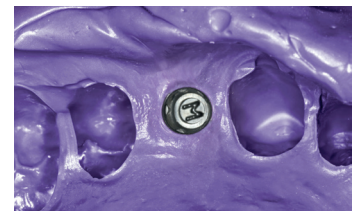


Polyether meets CAD/CAM workflow.





Impregum Super Quick Material is scannable with the current generation of scanners for a more efficient workflow.

3M™ Impregum™ Super Quick Polyether Impression Material

- First superfast polyether with 45 sec. working time and quick 2-minute setting
- Available in heavy, medium and light body consistencies
- For monophasic and 1-step impression techniques
- Perfect for impressions for single-unit crowns, implants, inlays/onlays and smaller bridges up to 3 units



The fastest way to polyether precision for smaller cases.

		Dispensing system	Impression technique*	Single-unit crown	Small implant (≤ 2 units)	Inlay/onlay (≤ 2 units)	Small bridge (≤ 3 units)	Recommended max. working time (min:sec)	Intra-oral setting time at 37°C (min:sec)
3M™ Impregum™ Penta™ Super Quick (Medium Body)		Monophase or 1-step	✓	✓	✓	✓	0:45	2:00	
3M™ Impregum™ Penta™ Super Quick (Heavy Body)		1-step	✓	✓	✓	✓	0:45	2:00	
3M™ Impregum™ Super Quick (Light Body)		1-step	✓	✓	✓	✓	0:45	2:00	

***Monophase Impression Technique:** one viscosity (medium body) in the tray and can also be syringed
1-step Impression Technique: two viscosities – either medium or heavy body in the tray and light body cartridge or syringe

Ordering information

Item no. Product information

NEW Upgrade Kits	
69404	3M™ Impregum™ Penta™ Super Quick Medium Body 1 base paste MB – 300 mL, 1 catalyst MB – 60 mL, 10 3M™ Penta™ Mixing Tips – Red
69411	3M™ Impregum™ Penta™ Super Quick Heavy Body / Light Body 1 base paste MB – 300 mL, 1 catalyst MB – 60 mL, 1 cartridge LB – 50 mL, 10 3M™ Penta™ Mixing Tips – Red, 5 3M™ Intra-oral Syringes Purple, 5 3M™ Garant™ Mixing Tips purple, 5 3M™ Garant™ Intra-oral Tips white
NEW Introductory Kits for 3M™ Pentamix™ 3 Automatic Mixing Unit	
69414	3M™ Impregum™ Penta™ Super Quick Medium Body 1 base paste MB – 300 mL, 1 catalyst MB – 60 mL, 13M™ Penta™ Cartridge MB, 10 3M™ Penta™ Mixing Tips – Red, 1 3M™ Penta™ Elastomer Syringe, 1 Polyether Tray Adhesive – 17 mL, 3 3M™ Astringent Retraction Paste capsules, 6 3M™ Dual Arch Trays
69415	3M™ Impregum™ Penta™ Super Quick Medium Body/Light Body 1 base paste MB – 300 mL, 1 catalyst MB – 60 mL, 2 cartridges LB – 50 mL each, 13M™ Penta™ Cartridge MB, 10 3M™ Penta™ Mixing Tips – Red, 5 3M™ Intra-oral Syringes Purple, 5 3M™ Garant™ Mixing Tips purple, 5 3M™ Garant™ Intra-oral Tips white, 1 Polyether Tray Adhesive – 17 mL, 3 3M™ Astringent Retraction Paste capsules, 6 3M™ Dual Arch Trays
69407	3M™ Impregum™ Penta™ Super Quick Heavy Body/Light Body 1 base paste HB – 300 mL, 1 catalyst HB – 60 mL, 2 cartridges LB – 50 mL each, 13M™ Penta™ Cartridge HB, 10 3M™ Penta™ Mixing Tips – Red, 5 3M™ Intra-oral Syringes Purple, 5 3M™ Garant™ Mixing Tips purple, 5 3M™ Garant™ Intra-oral Tips white, 1 3M™ Polyether Tray Adhesive – 17 mL, 3 3M™ Astringent Retraction Paste capsules, 6 3M™ Dual Arch Trays
NEW Refills	
69385	3M™ Impregum™ Penta™ Super Quick Medium Body 2 base pastes – 300 mL each, 2 catalysts – 60 mL each
69406	3M™ Impregum™ Penta™ Super Quick Heavy Body 2 base pastes HB – 300 mL each, 2 catalysts HB – 60 mL each
69379	3M™ Impregum™ Super Quick Light Body 4 cartridges – 50 mL each, 5 3M™ Intra-oral Syringes Purple, 5 3M™ Garant™ Mixing Tips purple

Item no. Product information

Cartridges for 3M™ Pentamix™ Automatic Mixing Units	
69387	3M™ Impregum™ Penta™ Super Quick Medium Body Cartridge
69401	3M™ Impregum™ Penta™ Super Quick Heavy Body Cartridge
71545	3M™ Universal Penta™ Cartridge for Pentamix™ 3
76994	3M™ Universal Penta™ Cartridge for Pentamix™ 2
77944	3M™ Penta™ Cartridge for 3M™ Pentamix™ Lite
Impressioning Trays	
71615	3M™ Directed Flow Impression Tray, Size S, Upper Jaw, 10 Trays
71616	3M™ Directed Flow Impression Tray, Size M, Upper Jaw, 10 Trays
71617	3M™ Directed Flow Impression Tray, Size L, Upper Jaw, 10 Trays
71618	3M™ Directed Flow Impression Tray, Size S, Lower Jaw, 10 Trays
71619	3M™ Directed Flow Impression Tray, Size M, Lower Jaw, 10 Trays
71620	3M™ Directed Flow Impression, Size L, Lower Jaw, 10 Trays
71622	3M™ Dual Arch Trays Posterior, 48 Trays
71623	3M™ Dual Arch Trays Anterior, 35 Trays
Related Products & Accessories	
71507	3M™ Intra-Oral Syringe – Purple, 20 Syringes
71508	3M™ Intra-Oral Syringes – Purple, 50 Syringes
77580	Cartridge Dispenser 1:1/2:1
77949	3M™ Penta™ Mixing Tips – Red, 50 Tips
71451	3M™ Polyether Intra-Oral Tips – White, 50 Tips
71454	3M™ Polyether Mixing Tips – Purple, 50 Tips
71461N	3M™ Garant™ Intra-oral Tips – White, 50 Tips
71210	3M™ Elastomer Syringe Set 1 Syringe; 4 Tips; 4 Barrels; 1 Brush
30600	3M™ Polyether Tray Adhesive – 17 mL
69408	3M™ Polyether Contact Tray Adhesive – 17mL
56944	3M™ Astringent Retraction Paste Refill Pack, 25 Capsules
56945	3M™ Astringent Retraction Paste Value Pack, 100 Capsules



3M Australia Pty Limited
 Building A, 1 Rivett Road
 North Ryde, NSW 2113
 1300 363 454
3M.com.au/ImpregumSQ

3M New Zealand Limited
 94 Apollo Drive
 Rosedale, Auckland 0632
 0800 808 182
3M.co.nz/ImpregumSQ

Not all products are available in both Australia and New Zealand. Always follow the Instructions for Use (IFU) and refer to IFU for full indications, precautions and warnings. Claims supported by 3M data on file. Comparative claims are in relation to 3M Impregum Soft Quick Step.

3M, “3M Science. Applied to Life.”, Garant, Impregum, Penta and Pentamix are trademarks of 3M. © 3M 2018. All rights reserved.