


3M™ RelyX™ Temporary Cements
A smart choice for temporary restorations

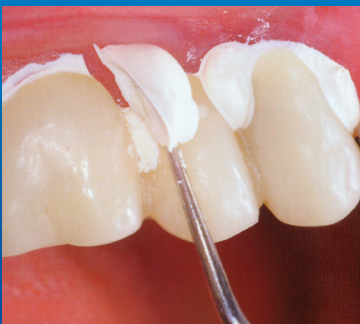

Cement is easily removed.

The ideal balance of reliable retention and easy removal

3M™ RelyX™ Temp NE Temporary Cement and 3M™ RelyX™ Temp E Temporary Cement offer you the best of both worlds—strong adhesion and easy removal for provisional cementation. Suitable for temporary cementation of provisionals and permanent restorations on implant abutments.

3M™ RelyX™ Temp NE Temporary Cement

A eugenol-free temporary cement that offers the ideal balance of reliable retention and easy removal for provisional restorations.

Indications

- Cementing temporary restorations
- Provisional cementation of permanent restorations
- Cementation of crown and bridges onto implant abutments

Product Benefits

- Extremely low film thickness helps ensure an optimal fit
- Strong adhesion offers high retention to the tooth, yet it is removed easily
- Excess cement can be easily removed
- Eugenol-free for universal indications because it will not inhibit the polymerization of composite resin cements; it's compatible with temporary crown and bridge materials, composite resin cements and composite resin core build-up materials

3M™ RelyX™ Temp E Temporary Cement

A temporary cement with a zinc-oxide, eugenol base for crown and bridge procedures.

Indications

- Cementing temporary restorations
- Provisional cementation of permanent restorations
- Cementation of crown and bridges onto implant abutments

Product Benefits

- Extremely low film thickness helps ensure an optimal fit
- Strong adhesion offers high retention to the tooth, yet it is removed easily
- Excess cement can be easily removed
- Has a sedative effect on the pulp due to the eugenol content



ITEM NO.	DESCRIPTION
56660	3M™ RelyX™ Temp NE Temporary Cement (non-eugenol) 1 – 30g Base Paste; 1 – 13g Catalyst Paste; 1 Mixing Pad (70144)



ITEM NO.	DESCRIPTION
35014	3M™ RelyX™ Temp E Temporary Cement (eugenol) 1 – 25g Base Paste; 1 – 18g Catalyst Paste; 1 Mixing Pad (70144)



3M Oral Care
2510 Conway Avenue
St. Paul, MN 55144-1000 USA

Phone 1-800-634-2249
Web 3M.com/dentalcements

3M Canada
Post Office Box 5757
London, Ontario N6A 4T1
Canada

Phone 1-888-363-3685

3M and RelyX are trademarks of 3M or 3M Deutschland GmbH. Used under license in Canada. All other trademarks are owned by their respective owners. Please recycle. © 3M 2019. All rights reserved.
70-2009-3461-3