



Clarity™

Aesthetic Orthodontic Solutions



3M™ Clarity™
ADVANCED Clear Braces

Brilliant smiles.
Exceptional performance.

Performance and aesthetics that exceed expectations.

3M™ Clarity™ Advanced Ceramic Brackets are a leap forward in orthodontic ceramic brackets. Through the combination of progressive technology and modern design, Clarity Advanced Ceramic Brackets deliver brilliant aesthetics, small design, predictable debonding and enhanced patient comfort. And, they're available with renowned APC™ Flash-Free Adhesive, to make the bonding experience faster and more efficient than ever before. These advanced features enable you to offer the ultimate combination of aesthetics and performance to your patients.

Attractive appearance. Distinct design.

Brilliant aesthetics.

Let your patient's smile shine through in every way with Clarity Advanced Ceramic Brackets. These ultra-small translucent brackets are designed to blend with one's natural tooth colour, while resisting staining and discolouration throughout treatment.




Trusted strength and small bracket size.


Showcasing a fine-grained ceramic material and created through an injection-molding process, Clarity™ Advanced Ceramic Brackets offer trusted strength in a small size and one of the lowest profiles in the industry. Additionally, Clarity Advanced Ceramic Brackets offer ample under tie-wing space, allowing for flexible treatment and ligation options.

Small grain size provides trusted strength in a smaller bracket.

Equal magnification.
Source: 3M lab.

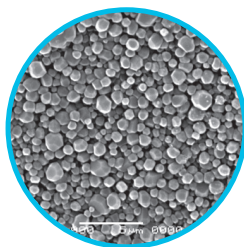

Fine-grained polycrystalline ceramic material
Clarity Advanced Ceramic Brackets

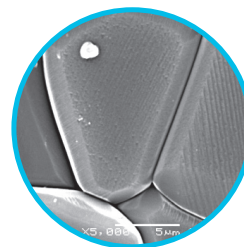

Coarse-grained monocrystalline ceramic material
3M Metal-reinforced brackets


Low anterior profile

0.0755 in.
(1.92 mm)



Clarity Advanced Ceramic Brackets
(Lower Anterior)

3M™ Clarity™ Advanced Ceramic Brackets: A revolution in ceramic bracket technology.

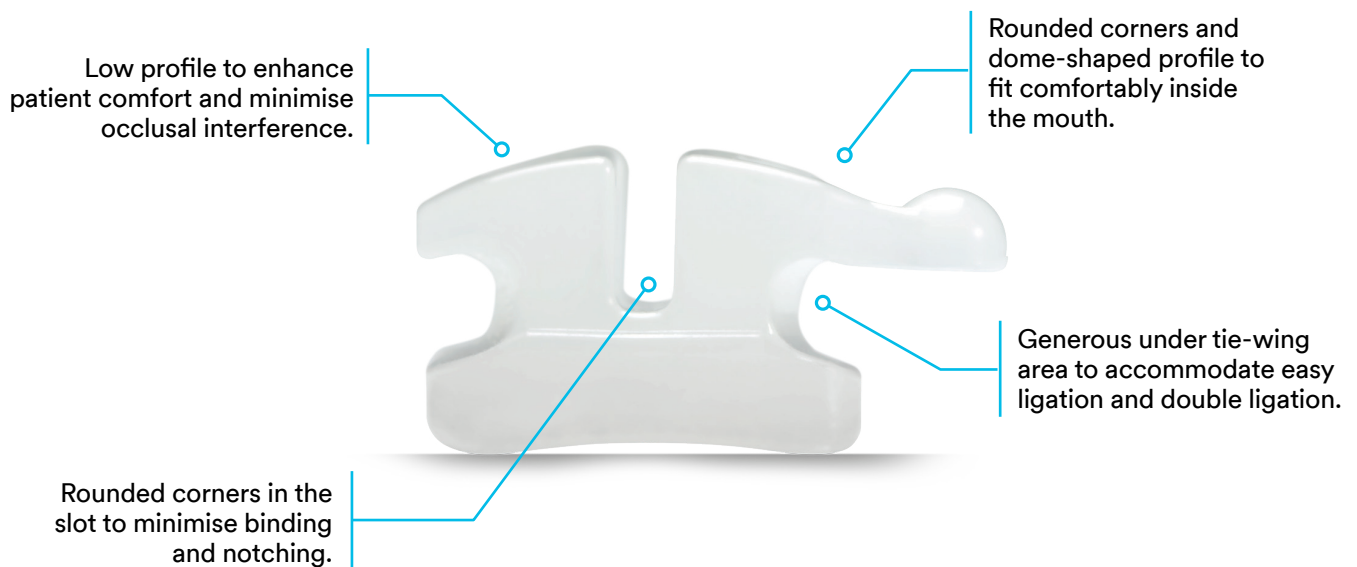

Treat aesthetically without compromise.



3M™ Clarity™ Advanced
Ceramic Brackets

Micro-crystalline, mechanical locking bonding surface for optimum bond strength.

Tooth-specific base geometry and compound contour base design for precise tooth fit.

Proprietary stress concentrator for predictable squeeze debonding.

Rounded bi-directional ball hook angled for patient comfort.

ID dimple with colour ink dots for bracket identification.

Base flange to facilitate easy bracket placement and adhesive flash clean-up.

Horizontal and vertical markers for visual reference during bracket placement.



3M™ Forsus™
Class II Correctors

Available in the
3M™ MBT™ Appliance
System, Roth* and High
Torque prescriptions.

*3M version of this prescription.
No endorsement by the Doctor
is implied.

Simple, efficient bonding process.

Powered by technology, the 3M™ APC™ Flash-Free Adhesive Coated Appliance System is one of the most efficient bonding systems in orthodontics. Our innovative design creates a unique adhesive bonding substrate on the base of the bracket that saves you time and reduces patient chair time.

Average bracket bond failure rate of less than 2%.*



Available precoated with 3M™ APC™ Flash-Free Adhesive for superior bonding confidence.



Saves Time and Cost

- ✓ Faster and more efficient bonding than non-precoated brackets.
- ✓ Eliminates the flash removal step from bonding.
- ✓ Simple, easy clean-up at debonding.



Enhanced Bonding Efficiency

- ✓ Bond time is reduced up to 40% per bracket.
- ✓ No flash removal reduces accidental bracket movement after positioning, reducing the need to reposition.
- ✓ Protects tooth enamel under adhesive.

*3M data on file.

Designed for happy patients and efficient practices.



Predictable, consistent debonding with proprietary stress concentrator.

Reduce the uncertainty of the ceramic bracket debonding process with our unique stress concentrator, a hallmark of the Clarity brand. This innovative design feature allows for mesial-distal debonding with a simple squeeze, so you can debond with confidence and help ensure the patients' debonding visit will be a positive part of the treatment experience.

Enhanced patient comfort with colours for any occasion.

Small, low profile bracket with smooth bracket surface and rounded corners enhance patient comfort. Patients can “blend in or be bold” by choosing from more than twenty clear or colourful 3M™ Alastik™ Easy-to-Tie Ligatures.





Clarity™

Aesthetic Orthodontic Solutions

Practice with clarity.

Our full portfolio of naturally-aesthetic treatment options are informed by 3M science and fueled by creativity. Innovation drives our leading-edge technology and delivers the choice and control you seek.

Our commitment is to collaborate with orthodontists to provide custom, flexible aesthetic treatment solutions that showcase your expertise and ensure precise outcomes. Together, we can change lives by creating healthy, beautiful smiles.



3M™ Clarity™ Aesthetic Orthodontic Solutions



3M™ Clarity™ Ultra Self-Ligating Brackets

Innovative, fully-aesthetic ceramic bracket designed for patient comfort and reliable results throughout treatment.



3M™ Clarity™ Advanced Ceramic Brackets

Clarity Advanced Brackets deliver brilliant aesthetics, small design, predictable debonding and enhanced patient comfort.

3M Oral Care

3M Australia Pty Limited
Building A, 1 Rivett Road
North Ryde, NSW 2113
Ph: 1300 363 484
3M.com.au/ortho

3M New Zealand Limited
94 Apollo Drive
Rosedale, Auckland 0632
Ph: 0800 441 622
3M.co.nz/ortho

3M, AlastiK, APC, Clarity, Forsus, MBT and Victory Series are trademarks of 3M.
© 3M 2019. All rights reserved.
70-2021-4916-0 1812