

Cementation Guide

Dentist:

Patient:

This is the restoration you are receiving and the respective pre-treatment (large box):

Z Zirconium Oxid

≤ 50 µm sandblast

Z1 Aluminium oxide 1 bar

Z2 Aluminium oxide 2 bar

ZR 3M™ Rocatec™

Clean, e.g. with air or alcohol

G Glass-Ceramic

Hydrofluoric 5-9,5%

G2 20 sec Hf

G6 60 sec Hf

G9 90 sec Hf

Clean, e.g. rinse with water and air-dry

C Composite

≤ 50 µm sandblast

Aluminium oxide 1 bar

Clean, e.g. with air or alcohol

M Metal

≤ 50 µm sandblast

Aluminium oxide 1 bar

Clean, e.g. with air or alcohol

Pre-treatment according to material recommendation already done by the dental technician: Yes No

<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
18	17	16	15	14	13	12	11	21	22	23	24	25	26	27	28
48	47	46	45	44	43	42	41	31	32	33	34	35	36	37	38
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

This is the cementation procedure we recommend (small box):

A Adhesive

3M™ RelyX™ Ultimate Cement



See step-by-step guide

SA Self-adhesive

3M™ RelyX™ Unicem (2) Cement



See step-by-step guide

C Conventional

3M™ Ketac™ Cem (Plus) Cement



See step-by-step guide


T Temporary

3M™ RelyX™ Temp NE Cement



Product information

Available resources

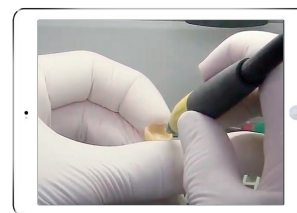

Step-by-step videos for cementation procedures

Please note that the final decision for the appropriate cementation procedure is the responsibility of the treating dental professional. Refer to relevant manufacturers Instructions for Use for complete information.

PLEASE CLEAN AFTER TRY IN ACCORDING TO MATERIAL RECOMMENDATION!

Note: _____ Date and signature: _____