3M Science. Applied to Life.™

Extended fluoride release varnish... because staying in contact matters.

Clinpro™ White Varnish

5% Sodium Fluoride
White Varnish with
Tri-Calcium Phosphate



Summary of advantages

- **Contains 22,600 ppm fluoride**
- **Adheres to teeth for continuous fluoride contact**
- **Provides a steady, sustained release of fluoride**
- **Migrates to "hard to reach" areas, including interproximal spaces**
- **Easy to apply to moist teeth**
- **Relieves hypersensitivity**
- **3 year shelf life**

Continuous fluoride contact.

Clinpro™ White Varnish 5% Sodium Fluoride with Tri-Calcium Phosphate (TCP) is a proprietary formula that contains fluoride, calcium and phosphate – components naturally found in saliva. It delivers targeted and sustained fluoride and calcium release.

Clinpro White Varnish flows across tooth surfaces and comes in contact with saliva, releasing fluoride and calcium over the life of the varnish – up to 24 hours!

Flows better.

Clinical studies show that Clinpro White Varnish flows to more than twice as many tooth surfaces after application.


Steady fluoride release.

Unlike varnishes that release most of their fluoride in the first hour, Clinpro White Varnish provides a prolonged, steady delivery of fluoride, allowing absorption into the tooth.


Independent studies have shown the value of sustained release of fluoride from varnishes:

"Most studies showing a higher efficacy of fluoride varnish attribute it to its increased contact time."*

*Source: Khattak MF, et al. Comparison of three topical fluorides using computer imaging. J Clin Pediatr Dent 2005; 30(2): 139-144.

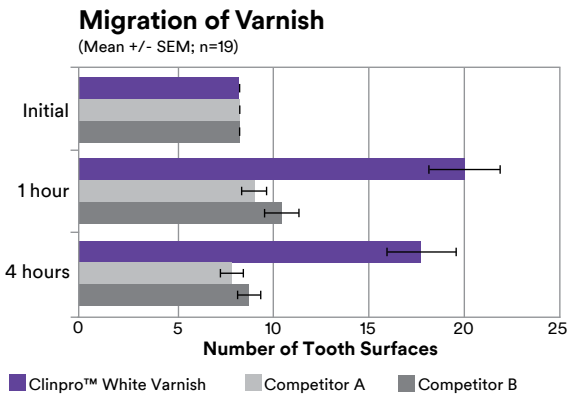

Clinpro™ White Varnish 5% Sodium Fluoride or the competitive varnishes were applied to 8 anterior maxillary surfaces in vivo. Migration of the varnishes was measured at 1 hour and 4 hours. Source: 3M internal data on file

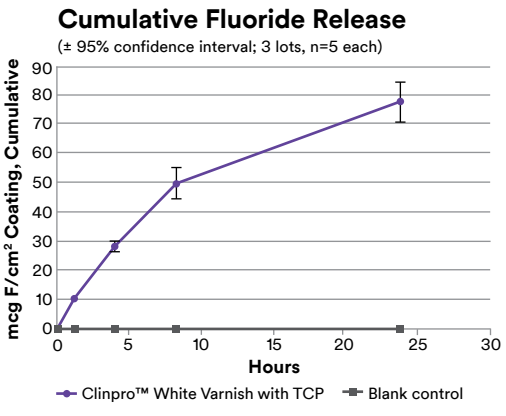


Clinpro™ White Varnish 5% Sodium Fluoride delivers a steady release of fluoride throughout the wear time. Source: 3M internal data on file

Easy to apply.

Simply apply one thin coat, sweeping the brush across tooth surfaces – and the Clinpro White Varnish 5% Sodium Fluoride formula does the rest.

- Sets in the presence of saliva – no drying needed.
- Enhanced flow characteristics allows Clinpro White Varnish to reach interproximal areas other varnishes may miss.
- Covers and adheres to the teeth until brushed off.
- Patient can eat and drink following application.

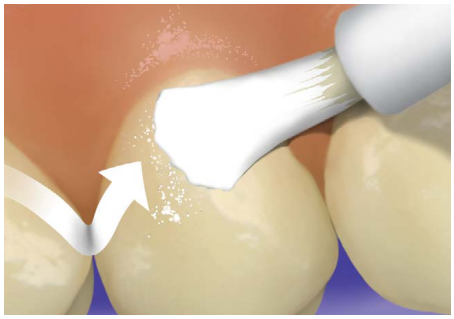

Apply Clinpro White Varnish 5% Sodium Fluoride evenly in a thin layer over treatment area(s) with sweeping, horizontal brush strokes. No suction required.

Protection you can't see.

Clinpro White Varnish is virtually invisible on the teeth – but leaves a temporary thin coating so patients know it's working.



Immediately after Clinpro White Varnish 5% Sodium Fluoride was applied.





Ordering Information

Item #	Product Information
	3M™ Clinpro™ 5% Sodium Fluoride White Varnish – 50 Pack Each 50 Pack includes: 50 unit-dose packets (0.5ml each) of one flavour and applicator brushes in individual foil pouches.
12249C	**Cherry Flavour**
12249L	**Melon Flavour**
12249	**Mint Flavour**
	3M™ Clinpro™ 5% Sodium Fluoride White Varnish – 100 Pack Each 100 Pack includes: 100 unit-dose packets (0.5ml each) of one flavour and 100 applicator brushes.
12250C	**Cherry Flavour**
12250L	**Melon Flavour**
12250	**Mint Flavour**

3M.com.au/PreventiveCare

Scan to find out more

3M Oral Care

3M Australia Pty Limited
Building A, 1 Rivett Road
North Ryde, NSW 2113
Ph: 1300 363 454
3M.com.au/oralcare

3M New Zealand Limited
94 Apollo Drive
Rosedale, Auckland 0632
Ph: 0800 80 81 82
3M.co.nz/oralcare

New Zealand, only Mint flavour available. Always follow the Instructions For Use. Claims supported by data on file. 3M Science. Applied to Life, 3M and Clinpro are trademarks of 3M. Please recycle. Printed in Australia.
© 3M 2021. All rights reserved. AV011477583.

FLAME_3M5970 12/20