



Littmann®
Stethoscopes

A new level of learning.


Quality makes a difference when you're learning auscultation skills to help you diagnose patients. Start with a 3M™ Littmann® Stethoscope to hear the subtle sounds you don't want to miss. It's built to last and helps you perform your best.



The outstanding acoustics of 3M™ Littmann® Cardiology and Digital Stethoscopes can help you diagnose patients and make decisions confidently—even with hard-to-hear sounds like aortic stenosis murmurs, S3 gallop rhythms and faint pulmonary anomalies.

Designed to provide excellent sound quality, these stethoscopes feature tunable diaphragm technology, so you can hear low and high frequency sounds without repositioning or turning over the chestpiece. The dual-lumen tubing eliminates the rubbing noise of traditional double tubes, and the adjustable headsets and snap-tight soft-sealing eartips ensure an excellent acoustic seal and comfortable fit.

The 3M™ Littmann® CORE Digital Stethoscope takes learning to a new level. It connects to Eko software, allowing you to visualize heart sound waveforms in real time.

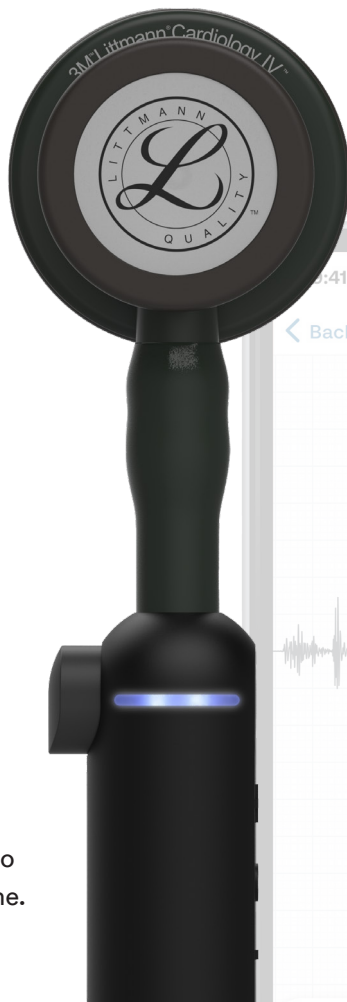

3M™ Littmann® CORE Digital Stethoscope

Up to 40x amplification (at peak frequency, vs. analog mode) and in-app sound wave visualization for a new level of learning.



3M™ Littmann® Cardiology IV™ Stethoscope

Outstanding acoustics for detailed diagnostic auscultation with the versatility of adult and pediatric sides of the chestpiece.

♥ 67 BPM

3M™ Littmann® Cardiology IV™ Stethoscope



Standard stainless steel



6151 (22") Black 6152 Black 6154 Navy 6155 Hunter 6156 Plum



6158 Raspberry 6159 Rose Pink 6184 Burgundy

Black/black



6163 Black 6168 Navy 6200 Black 6201 Black 6203 Black

Mirror



6177 Black 6170 Burgundy

Rainbow/black



6165 Black 6205 Plum

High-polish smoke/smoke



6238 Gray

High-polish rainbow/smoke



6240 Black 6241 Raspberry

High-polish rainbow/black



6242 Navy

High polish smoke/black



6202 Navy 6204 Black 6232 Black 6234 Caribbean

Champagne/black



6176 Burgundy 6190 Caribbean

Champagne/smoke





6179 Black

Product Features	
Acoustic Performance	9 (1-10 scale)
Chestpiece Design/Technology	Machined stainless steel, two-sided
Headset Material	Lightweight aluminum alloy
Length/Total Weight	22"/56 cm/167 grams, 27"/69 cm/177 grams
Diaphragm Details/Diameter	Two, single-piece tunable diaphragms 1.7"/4,3 cm, 1.4"/3,4 cm
Warranty	7 years



3M Health Care
3M Medical Solutions Division
2510 Conway Ave., St. Paul, MN 55144 USA
1-800-228-3957 | Littmann.com

3M™ Littmann® CORE Digital Stethoscope

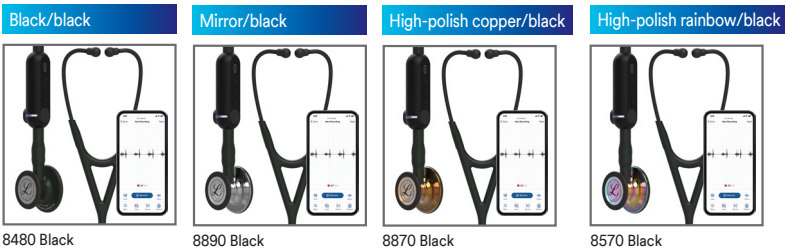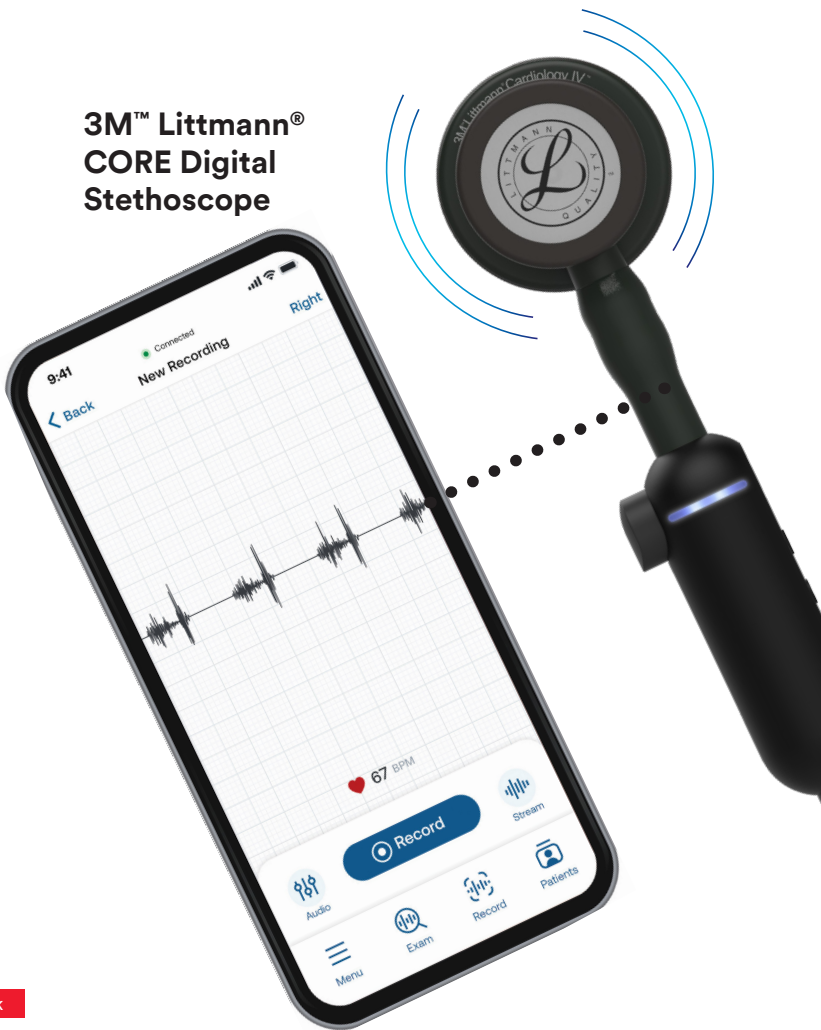

8480 Black 8890 Black 8870 Black 8570 Black

Product Features	
Acoustic Performance	10 (1-10 scale)
Chestpiece Design/Technology	Machined stainless steel, two-sided
Headset Material	Lightweight aluminum alloy
Length/Total Weight	27"/69cm/232g
Diaphragm Details/Diameter	Two, single-piece tunable diaphragms 1.7"/4,3 cm, 1.4"/3,4 cm
Warranty	2 years

*Colors shown are approximate.

Explore Littmann stethoscopes at [Littmann.com](https://littmann.com) or call the 3M Health Care Helpline at 1-800-228-3957.

3M, Littmann, the L Littmann logo, and Cardiology IV are trademarks of 3M. Eko trademarks, including CORE, are marks and/or registered marks of Eko Devices, Inc. Used under license in Canada. Please recycle. Printed in U.S.A. © 3M 2023. 70-2011-8178-4