

Adper™ Easy Bond

Self-Etch Adhesive



ONE bottle. ONE coat.
One easy, strong bond.

Summary of Advantages

- One-bottle, one-coat convenience
- Bond performance equal to or greater than other one-bottle self-etch adhesives
- Reduced risk of post-operative sensitivity
- Faster procedure for less risk of contamination
- Moisture tolerant – works on wet or dry enamel and dentin
- Low microleakage
- Available in unit dose delivery

Use it because it's easy. Use it because it's strong.

Unlike total-etch adhesive systems, Adper Easy Bond adhesive simultaneously etches and bonds to dentin and enamel in one step. This means:

- You only need one bottle, one coat and one step to achieve a strong, reliable bond.
- Fast, easy application with no mixing required.

Competitive Products Comparison

	Adper™ Easy Bond	Xeno® IV	iBond™
Higher bond strengths	✓		
One-bottle, one-coat convenience	✓	2 Coats	3 Coats
Ethanol-based solvent system – less volatile than acetone	✓	Acetone	Acetone
Vitrebond™ copolymer for more consistent performance under variable humidity	✓		
Bonded filler for improved bond strength	✓		

Source: 3M ESPE internal data




Performance backed by 25 years of adhesive experience.

Adper Easy Bond adhesive builds on much of the same chemistry as Adper™ Scotchbond™ Multi-Purpose Adhesive and Adper™ Single Bond Plus Adhesive. So while you may choose this one-coat self-etch for convenience, you can use it with confidence knowing it's built upon a strong tradition of bonding technology and performance.

Adper™ Easy Bond Self-Etch Adhesive offers high bond strength to cut enamel and dentin.

Adper Easy Bond adhesive utilizes a filler technology for optimized bond performance. Independent studies prove its bond strength is equal to or greater than that of one-bottle, and compares favorably to two-bottle self-etch adhesives.

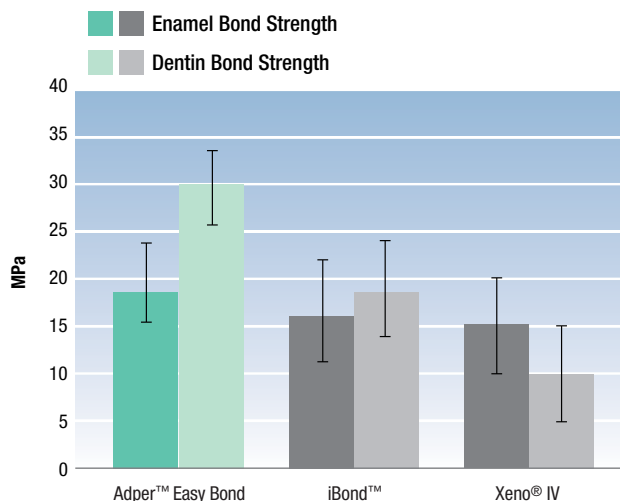

Bond strength to cut enamel and dentin.

Source: 3M ESPE internal data

One coat means a faster application procedure.


The following illustrates how simple it can be using Adper Easy Bond adhesive for direct light-curing restoration, whether using our convenient, disposable Unit Dose applicator or bottle.



Step 1. Apply adhesive to tooth surface for a total of 20 seconds.



Step 2. Dry the adhesive for 5 seconds.

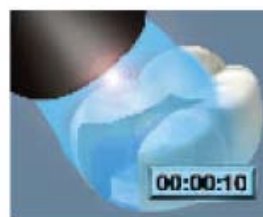

Step 3. Light cure for 10 seconds.

Ordering Information



Item #	Product Information
41245	Adper™ Easy Bond Self-Etch Adhesive Intro Kit Kit includes: 5ml adhesive vial; 50 disposable mixing wells; 50 disposable microbrushes
41246	Adper™ Easy Bond Self-Etch Adhesive Refill 5ml adhesive vial
41248	Adper™ Easy Bond Self-Etch Adhesive Unit Dose Intro Pack 40 unit dose
41249	Adper™ Easy Bond Self-Etch Adhesive Unit Dose Giant Pack 100 unit dose
41252	Adper™ Easy Bond Self-Etch Adhesive Unit Dose Bulk Pack 200 unit dose

Customer Care Center 1-800-634-2249

www.3MESPE.com

3M ESPE

Dental Products

3M Center
Building 275-2SE-03
St. Paul, MN 55144-1000
U.S.A.
1-800-634-2249

3M Canada

Post Office Box 5757
London, Ontario N6A 4T1
Canada
1-888-363-3685

Please recycle. Printed in U.S.A.
© 3M 2007. All rights reserved.
70-2009-3961-2

3M, ESPE, Adper, Scotchbond and Vitrebond are trademarks of 3M or 3M ESPE AG. Used under license in Canada. iBond and Xeno are not trademarks of 3M.