

3M


Tegaderm™

Chlorhexidine Gluconate (CHG)
I.V. Securement Dressing

The difference is clear.

Only Tegaderm™ CHG Dressings combine the most critical elements of I.V. site care into a single, easy-to-use product.

- Full site visibility
- Antimicrobial CHG protection
- Catheter securement
- Breathability

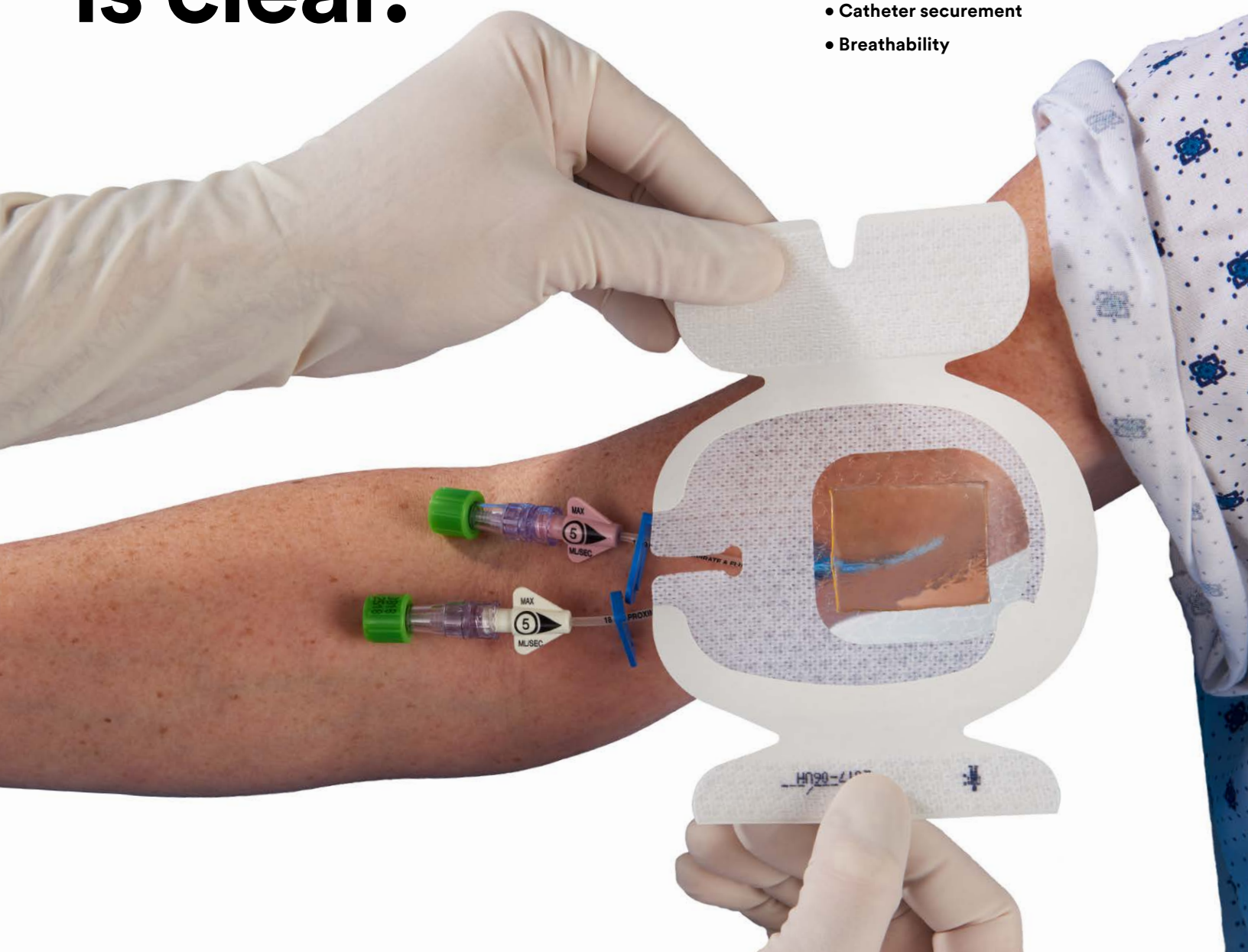

1657, 1658, 1659

3M™ Tegaderm™ Chlorhexidine Gluconate (CHG) I.V. Securement Dressing


1657, 1658, 1659

Integrated antimicrobial (CHG) gel pad

- 2% CHG is integrated into the dressing's gel pad, so it can't be dropped, forgotten or applied incorrectly
- Transparent gel pad provides full site visibility
- Conforms around catheter and hub for complete CHG coverage
- Delivers immediate and continuous CHG protection for up to 7 days, even in the presence of blood and exudate
- Absorbs blood and other fluids
- Activates without requiring skin moisture

Large securement tape strip with notch

- Promotes consistent application and enhances stabilization
- Includes adhesive-free tabs to minimize potential to stick to gloves or to itself

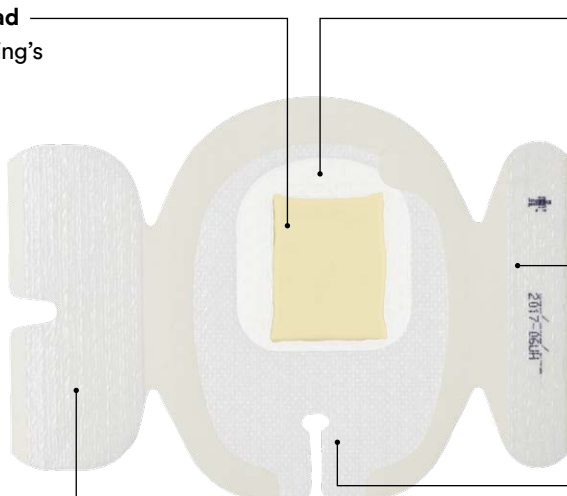

Highly breathable transparent film

- Provides a barrier to external contaminants including liquids, bacteria and viruses*
- Promotes moisture evaporation

Reinforcing documentation tape strip

- Includes pre-printed labels for documenting dressing changes
- Film coating provides a barrier to external contaminants*
- Can also be used to secure hubs, lumens and tubing

Film-coated stabilization border

- Keyhole notch conforms around catheter hubs and lumens
- Conforming edge border flexes with patient movement and may help to reduce edge lift
- Patterned-coated film adhesive holds strongly, manages moisture and releases gently

PICC



Jugular



Subclavian






Arterial



Femoral



Product	Product Number	CHG Gel Pad Size		Overall Dressing Size		Units/Box	Boxes/Case
3M™ Tegaderm™ Chlorhexidine Gluconate (CHG) I.V. Securement Dressing							
	1657	1½ in x 1¾ in	3 cm x 4 cm	3½ in x 4½ in	8,5 cm x 11,5 cm	25	4
	1658	1½ in x 1¾ in	3 cm x 4 cm	4 in x 4¾ in	10 cm x 12 cm	25	4
	1659	1½ in x 2¼ in	3 cm x 7 cm	4 in x 6½ in	10 cm x 15,5 cm	25	4

To learn more about Tegaderm™ CHG Dressings or to request a product evaluation at your facility, visit us at 3M.com/TegadermCHG, contact your 3M Critical & Chronic Care Solutions representative or call the 3M Health Care Customer Helpline at 1-800-228-3957. Outside of the United States, contact the local 3M subsidiary.

* *In vitro* testing shows that transparent film provides a viral barrier from viruses 27nm in diameter (e.g. HCV) or larger (e.g. HBV and HIV) while the dressing remains intact without leakage.



Critical & Chronic Care Solutions Division

3M Health Care
2510 Conway Avenue
St. Paul, MN 55144 USA

Phone 1-800-228-3957
Web 3M.com/TegadermCHG

3M Canada

P.O. Box 5757
London, Ontario N6A 4T1
Canada

Phone 1-800-364-3577
Web 3M.ca/HealthCare

Please recycle. Printed in U.S.A.
© 3M 2011, 2016. All rights reserved.
3M and Tegaderm are trademarks of 3M.
70-2009-9594-5