

3M™ Attest™ Rapid 5 Steam-Plus Test Pack

Includes 3M™ Attest™ Rapid Readout Biological Indicator
3M™ Comply™ (SteriGage™) Steam Chemical Integrator (AAMI Class 5)

#1 Utilized Steam Sterilization Test Pack
in U.S. Hospitals*

Fast, Easy to Use Steam Packs

Comply with Standards and Reduce Worry

- 3M™ Attest™ Rapid 5 Steam-Plus Test Pack is a process challenge device (PCD) that contains a 3-hour 3M™ Attest™ Rapid Readout Biological Indicator plus a 3M™ Comply™ (SteriGage™) Steam Chemical Integrator (AAMI Class 5).
- Compliant with AAMI and AORN standards that recommend a BI PCD containing a Class 5 Integrating Indicator be used in every load containing implants.

Fast Turnaround Times for Implant Load Releases

- 3-hour rapid readout BI provides an 88% decrease in waiting times for biological results compared to 24 hour BIs.

Helps Reduce Risks

- With faster results, decisions can be made and action taken prior to surgery, reducing the risks and cost associated with surgical infections.

Helps Reduce Recalls

- When each load is monitored and quarantined until BI results are available, the severity of recalls may be reduced.



*2010 HPIS Industry Data

3M Infection Prevention Solutions

Innovation
On A Mission

3M

Easy to Use Process Challenge Device and Monitoring System




3M™ Attest™ 1292 Rapid Readout Biological Indicator

- Contains a large number of highly resistant spores
- Provides the fast readout times (3-hour)
- BIs are the only sterilization process monitoring device that provide a direct measure of lethality



3M™ Attest™ 290 Auto-reader for Steam

- BI results obtained using the 3M™ Attest™ 290 Auto-reader
- Does not require time consuming calibration
- Automated results

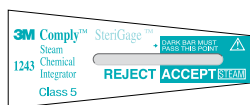

3M™ Comply™ (SteriGage™) Steam Chemical Integrator (AAMI Class 5)

- Conforms to ANSI/AAMI/ISO 11140-1:2005 Class 5 Integrating Indicator performance requirements
- Referenced in AAMI ST79 recommended practices for implant load release in emergency situations only
- Does not require color interpretation



Every Load Monitoring Helps Reduce Recalls

- When every load is monitored with a 3M Attest Rapid 5 Steam-Plus Test Pack, the severity of recalls may be reduced.
- Join the hundreds of other healthcare facilities in the United States utilizing the 3M Attest Rapid 5 Steam-Plus Test Pack to monitor every load and ensure that every patient receives a consistent, high level of care.



Ordering Information

Catalog No.	Description	Packaging	Shipping
41382	3M Attest Rapid 5 Steam-Plus Test Pack – monitors 132°C (270°F), ≥4 minute vacuum-assisted and 121°C (250°F), ≥30 minute gravity displacement steam sterilization processes (Must be used with 3M Attest 290 Auto-reader)	16 each/box 16 controls/box	4 boxes/case
41382F	3M Attest Rapid 5 Steam-Plus Test Pack	16 each/box 4 controls/box	4 boxes/case
290	3M Attest Auto-reader for Steam, for use with 3M Attest Rapid Readout Biological Indicators for Steam 1291 and 1292	1 unit/box	1 box/case

For more information contact your 3M Sales Representative, call the 3M Health Care Customer Helpline at 1-800-228-3957 or visit www.3M.com/infectionprevention



Infection Prevention Division
3M Health Care
3M Center, Building 275-4E-01
St. Paul, MN 55144-1000
U.S.A.
1 800 228-3957
www.3M.com/infectionprevention

Please recycle. Printed in U.S.A.
3M, Attest, Comply and SteriGage are trademarks of 3M.
© 3M 2011. All rights reserved.
Revised 4/2011
70-2009-7453-6 (411) DPI