


3M™ Bair Hugger™ Therapy
Lithotomy Underbody Blanket
Model 585

ACHIEVE
SCIP-Inf-10
GOALS



Warm From Below

The lithotomy underbody blanket is positioned on the procedure table before the patient arrives at the room. Designed to deliver full, unrestricted patient access and flexibility for procedures involving the lower extremities and abdominal, peritoneal and pelvic cavities.

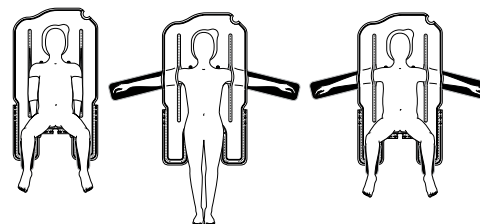

3M Infection Prevention Solutions

Innovation
On A Mission

3M

Features

- Underbody design provides full, unrestricted access to the patient and accommodates lateral, supine and lithotomy positions
- Fluid outlets minimize pooling of fluids on the surface of the blanket
- Pass-through slit design allows for flexible patient positioning and accommodates the use of a drawsheet
- Generously-sized removable peritoneal cut-out at the base of the blanket provides clinician unobstructed access to patient
- Integrated tie strips can be loosely secured around the stirrups/leg supports
- The blanket draws up next to the patient when inflated to surround them in warmth
- Tape and tuck flaps secure the blanket to the O.R. table


Ordering Information

Product	Part Number	Specifications
Lithotomy Underbody Blanket Model 585	58501	Size: 64" x 36" (163 x 91 cm) Drape: 24" x 24" (61 x 61 cm) Weight: 5 oz (153 g) Units/Case: 10


To order, please contact your local 3M Patient Warming sales representative or call **Customer Service at 1-800-733-7775**.

Product Offerings


Preoperative Blankets

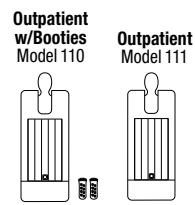


Postoperative Blankets

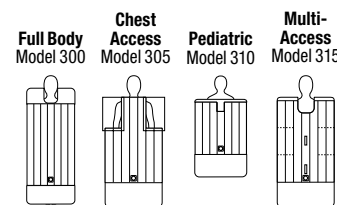


Specialty and Cardiac Blankets

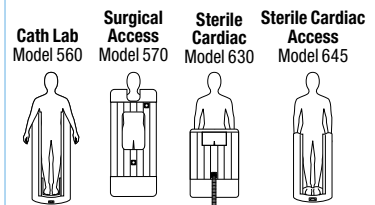


Intraoperative Blankets

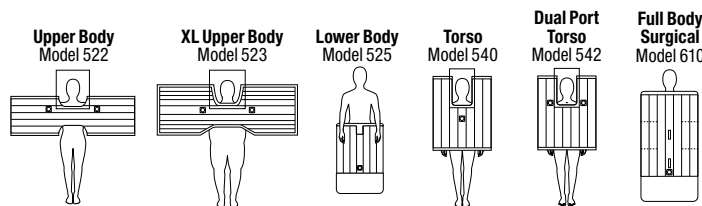

Pediatric Blankets

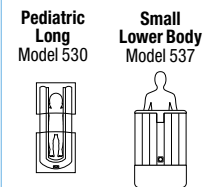

Underbody Series Blankets

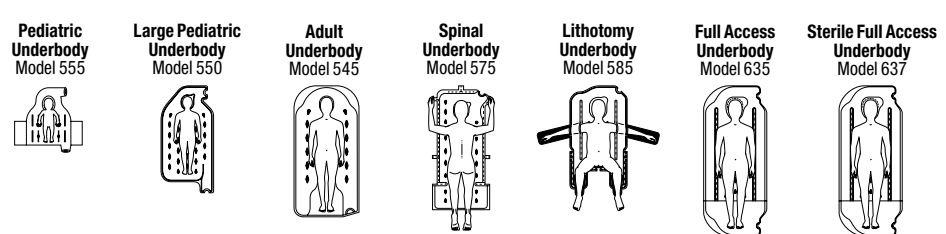

Arizant Healthcare Inc.,
a 3M company
10393 West 70th St.
Eden Prairie, MN 55344 USA
Phone 800-733-7775
Fax 800-775-0002
www.bairhugger.com

3M is a trademark of 3M Company, used under license in Canada. BAIR HUGGER and the BAIR HUGGER logo are trademarks of Arizant Healthcare Inc., used under license in Canada. ©2012 Arizant Healthcare Inc. All rights reserved. 601907F 7/12